


Simple, Reliable and Fast Electrolytic Detachment with Real-time Feedback

- **Place** NumenFR™ Detachment System over Numen™ coil's proximal end of delivery wire, until the **System Ready Indicator** illuminates green and one short beep sounds.
- **Press** the **Detachment Button** to activate the detachment.
- **Observe** the **Status Indicator** – a flashing green light indicates detachment in process – for 5 -10 seconds until it illuminates solid green light with 3 short beeps or 1 long beep sound.
- **Confirm** successful detachment under fluoroscopy and pull out NumenFR™ Detachment System.

Description	Product Catalog Number
NumenFR™ Detachment System	E101

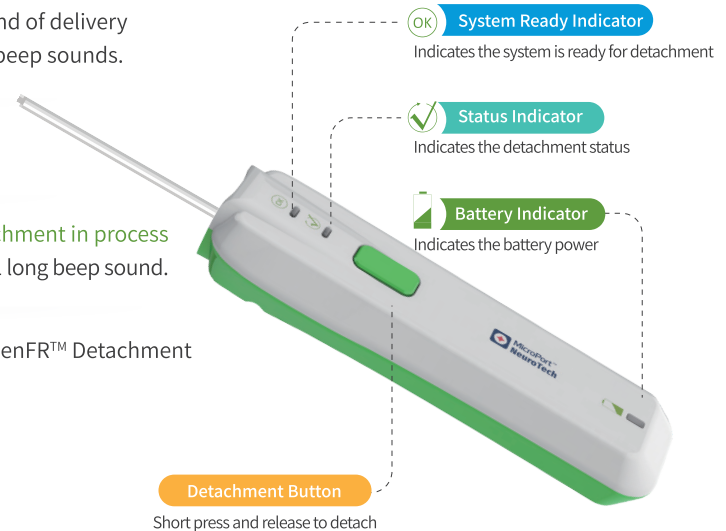

Caution: Do not use detachment systems other than NumenFR™ Detachment System

NumenFR™ Indicator Modes and Suggested Actions



INDICATION	MODE	SUGGESTED ACTION
OK SOLID GREEN 1 SHORT BEEP	Ready for detachment	Press the Detachment Button
✓ FLASHING GREEN	Detachment in progress	Wait for detachment
✓ SOLID GREEN 3 SHORT BEEPS/ 1 LONG BEEP	Detachment cycle complete	Confirm successful detachment under fluoroscopy. If the coil is not detached, repeat the detachment procedure
OK OFF	Connection error	Pull out NumenFR™ and place it over the delivery wire again
🔋 FLASHING YELLOW	Battery low	Replace the NumenFR™
OK ✓ 🔋 FLASHING YELLOW CONTINUOUS SHORT BEEPS	Battery error	
OK ✓ 🔋 FLASHING YELLOW CONTINUOUS SHORT BEEPS	System error	
OK 🔋 FLASHING YELLOW CONTINUOUS SHORT BEEPS	Short circuit of delivery wire	Replace the Numen™ Coil