

Superior Results with 3T and NavPro Technology

Transition-Free Design

Continuous tip movement reduces vascular trauma

0.41 mm Tapered Profile


Effortless entry into complex tortuous anatomy

Tri-Layered Atraumatic Tip

Flexible, durable, and safe for optimal support and stability

NavPro Technology

Smooth navigation and durable performance, with enhanced pushability and kink resistance

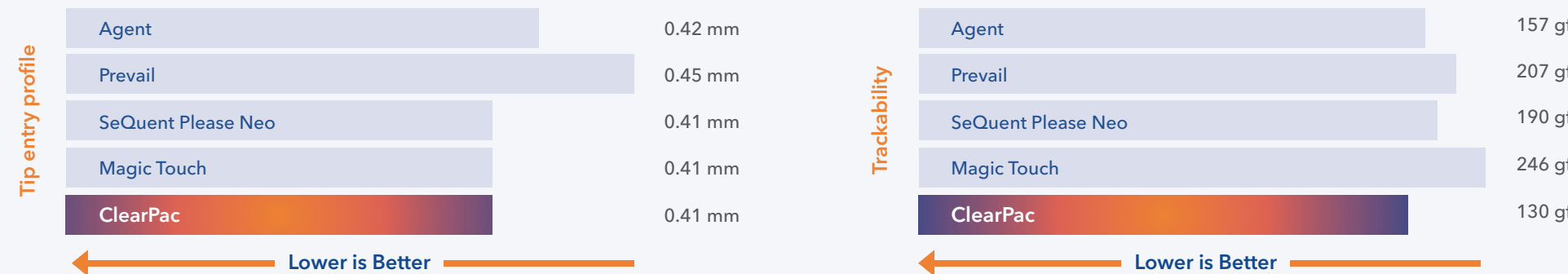

** Data on file; gf - Gram Force; mm - Millimeter
Agent is a registered brand of Boston Scientific. Prevail is a registered brand of Medtronic. SeQuent Please Neo is a registered brand of B. Braun. Magic Touch is a registered brand of Concept Medical.

Product Overview

Ordering Matrix

Balloon Length (mm)	Balloon Diameter (mm)							
	2.00	2.25	2.50	2.75	3.00	3.50	4.00	4.50
08	CLP20008	CLP22508	CLP25008	CLP27508	CLP30008	CLP35008	CLP40008	CLP45008
12	CLP20012	CLP22512	CLP25012	CLP27512	CLP30012	CLP35012	CLP40012	CLP45012
16	CLP20016	CLP22516	CLP25016	CLP27516	CLP30016	CLP35016	CLP40016	CLP45016
20	CLP20020	CLP22520	CLP25020	CLP27520	CLP30020	CLP35020	CLP40020	CLP45020
24	CLP20024	CLP22524	CLP25024	CLP27524	CLP30024	CLP35024	CLP40024	CLP45024
28	CLP20028	CLP22528	CLP25028	CLP27528	CLP30028	CLP35028	CLP40028	CLP45028
32	CLP20032	CLP22532	CLP25032	CLP27532	CLP30032	CLP35032	CLP40032	CLP45032
36	CLP20036	CLP22536	CLP25036	CLP27536	CLP30036	CLP35036	CLP40036	CLP45036
40	CLP20040	CLP22540	CLP25040	CLP27540	CLP30040	CLP35040	CLP40040	CLP45040
44	CLP20044	CLP22544	CLP25044	CLP27544	CLP30044	CLP35044	CLP40044	CLP45044
48	CLP20048	CLP22548	CLP25048	CLP27548	CLP30048	CLP35048	CLP40048	CLP45048
52	CLP20052	CLP22552	CLP25052	CLP27552	CLP30052	CLP35052	CLP40052	CLP45052

Technical Specifications

Available Balloon lengths (mm)	08, 12, 16, 20, 24, 28, 32, 36, 40, 44, 48, 52
Available Balloon Diameters (mm)	2.00, 2.25, 2.50, 2.75, 3.00, 3.50, 4.00, 4.50
Drug Dose	2.30 µg/mm ²
Balloon Material	Polyamide
Catheter Design	Rapid Exchange Delivery System
Usable Catheter Length (cm)	140
Guidewire Compatibility	0.014"

Fr - French; atm - Atmospheric Pressure

Guiding Catheter Compatibility	5 Fr (Ø2.00 – 4.00 mm) 6 Fr (Ø4.50 mm)
Catheter Shaft Diameter	1.95 Fr Proximal, 2.70 Fr Distal
Balloon Folds	3 Folds
Coating on Delivery System	Hydrophilic Coating
Balloon Nominal Pressure (NP)	8 atm
Rated Burst Pressure (RBP)	16 atm (Ø2.00 – 4.0 mm) 14 atm (Ø 4.50 mm)



Manufacturing Unit
Plot No. 105, Solitaire Industrial Park, Village Dehmikalan, Teh Sanganer, Jaipur, Rajasthan, India 303 007.
Manufacturing License No: MFG/MD/2022/000680

Head Office
Plot No. 121/F, Hosur Rd, 1st Phase, Bommasandra Industrial Area, Bengaluru, Karnataka 560099.

Marketing Office
Advant Navis Business Park Plot #7, Office No:1012, B Wing,10th Floor, Tower B, Noida-Greater Noida Expy, Sector 142, Noida, Uttar Pradesh, India 201 305.

International Office
Innvolution Healthcare Global Pte Ltd.
Marina Bay Financial Centre Tower 1, Level 11, 8 Marina Boulevard, Singapore 018 981.

Email: info@innvolution.com
Website: www.innvolution.com



Optimised Drug Delivery, Maximising Impact



ClearPac
Paclitaxel Coated Balloon

Powered by **ShellPac Coating Technology**

Optimised Drug Dose for Long-Term Arterial Healing



Paclitaxel

Lipophilic and anti-proliferative agent that prevents restenosis

Shellac - USFDA Approved

A hydrophilic resin enabling optimum release of Paclitaxel



Seamless 360° Shellac Coating
Prevents uneven drug loss while ensuring precise adhesion and release.




2.30 µg/mm² Paclitaxel effectively inhibits cell proliferation and prevents restenosis within the optimal dosage range.

Therapeutic Range of Paclitaxel* (2.00 - 3.50 µg/mm²)


BRAND	DOSE
Agent	2.00 µg/mm ²
Prevail	3.50 µg/mm ²
SeQuent Please Neo	3.50 µg/mm ²
ClearPac	2.30 µg/mm ²

Optimised Dose of Paclitaxel


Engineered for Efficient Drug Release with ShellPac Coating Technology

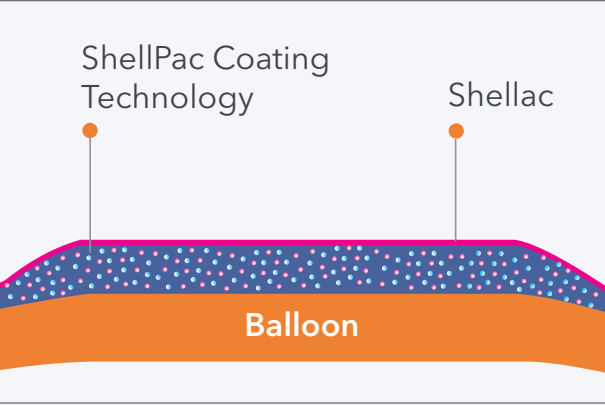


- ShellPac consists of a unique **1:1 formulation** of Paclitaxel and Shellac
- Shellac applied through a **360° spray coating** for uniform drug delivery
- Additional protective Shellac layer **enhance coating stability** and prevents early washout and drug loss during transit

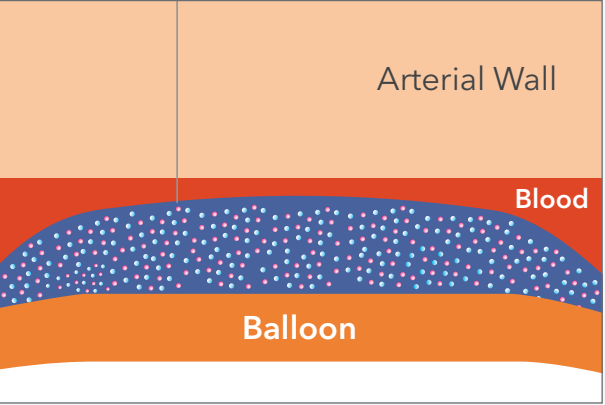

- **Hydrophilic Shellac** matrix swells upon contact with blood
- Controlled and optimised release of Paclitaxel during inflation

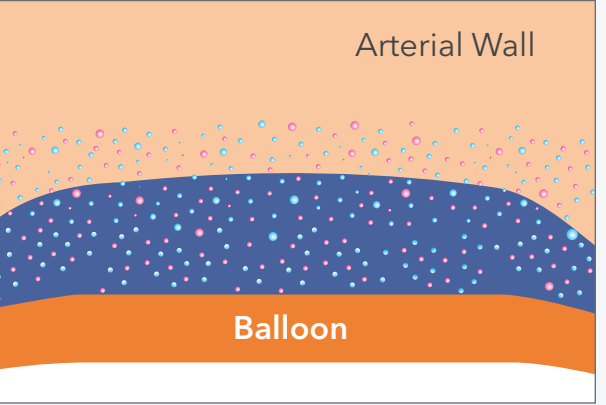

- Drug Release:**
- Quick and effective delivery of Paclitaxel within **45 secs**
 - Prevents excessive tissue growth
 - Reduces chances of **restenosis**

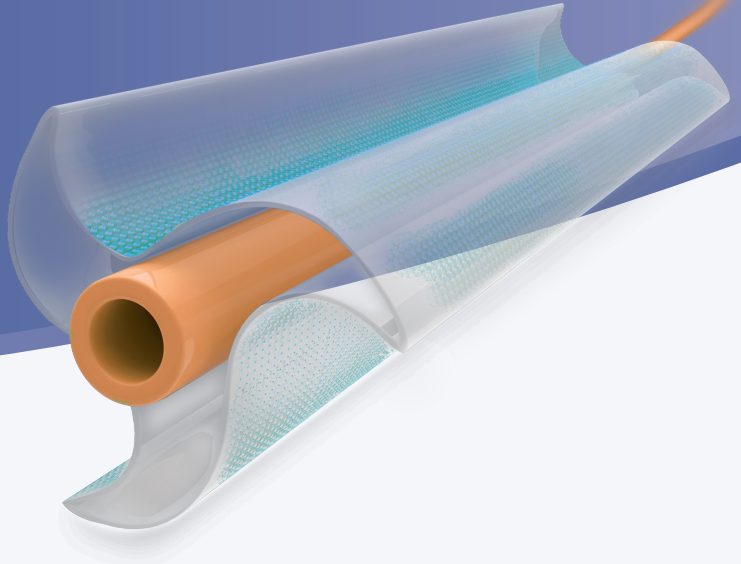

- **Efficient Drug Carrier**
Ensures uniform Paclitaxel distribution for targeted therapeutic action
- **Controlled Release**
Regulates Paclitaxel elution, maximising local efficacy
- **Strong Adhesion**
Prevents drug loss, ensuring effective delivery during inflation
- **Biocompatible & Safe**
A natural resin with excellent tolerance, reducing the risk of adverse reactions.
- **Enhanced Drug Stability**
Protects Paclitaxel from early washout, preserving its potency for sustained impact

Designed for:

- De-Novo Lesions
- Diffuse Lesions
- Bifurcation Lesions
- Ostial Lesions
- Calcified Lesions
- Long Lesions
- Small Vessels
- In-Stent Restenosis
- ST Elevation Myocardial Infarction
- Acute Coronary Syndrome
- High-Bleeding Risk
- Young Patients

*Ref - Future Pharmacol. 2024, 4, 775-787; **Data on file
Agent is a registered brand of Boston Scientific. Prevail is a registered brand of Medtronic. SeQuent Please Neo is a registered brand of B. Braun.