


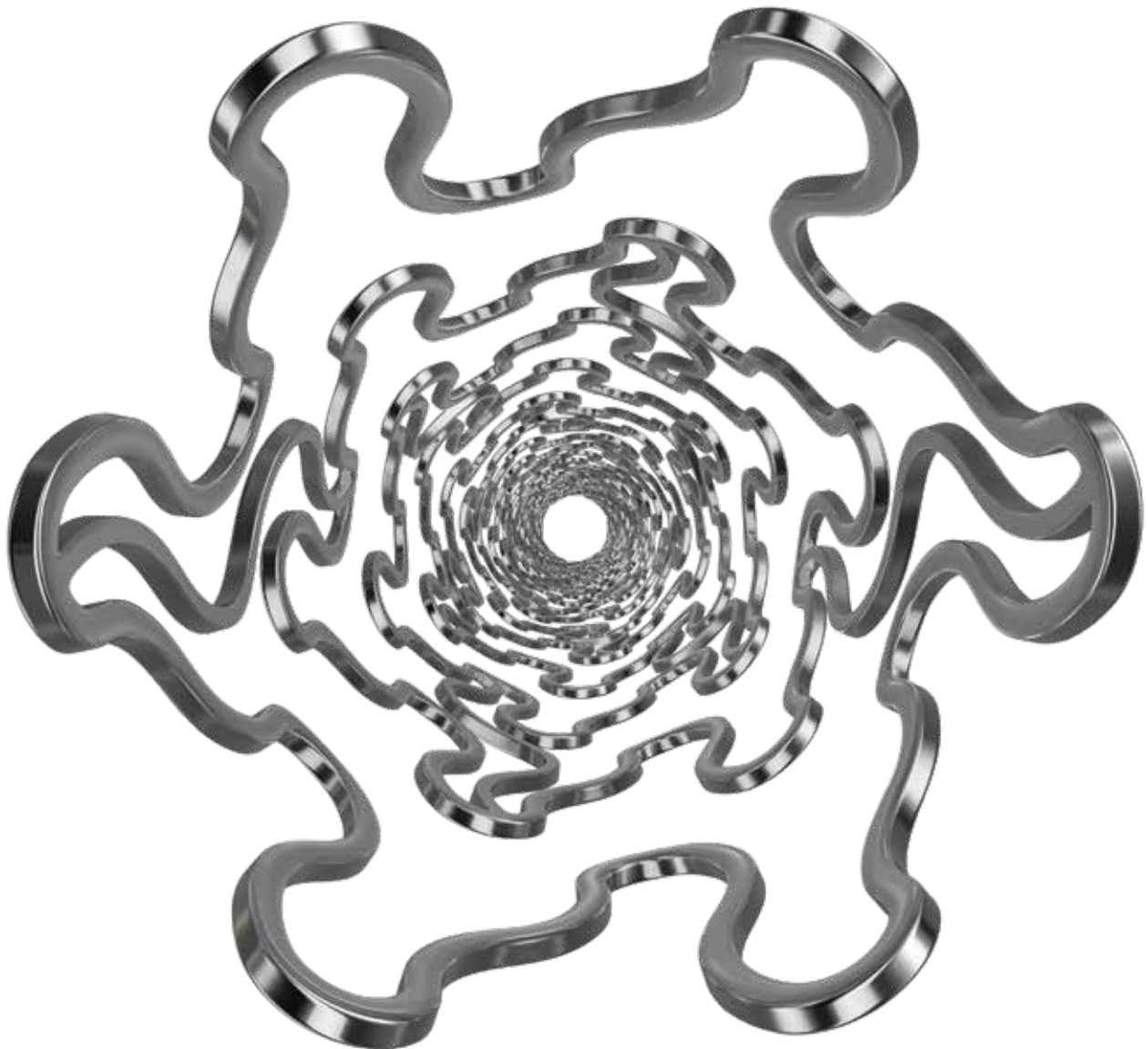

ETERNIA

**Everolimus Eluting
Coronary Stent System**

**Unblock with
Confidence!**

Restore Circulation with Eternia

What makes it the Future
of Coronary Stenting?

Everolimus

Enabling rapid
absorption into the
arterial wall

60 μ m Thin
Struts

1 mm Crossing
Profile

Improved deliverability
with minimal injury

L605 Cobalt
Chromium


High Radial Strength
(Force/Axial Length)
of 1.50 N/mm

Next-generation DES

Ensuring the well-being of your CAD patients!

Restore Circulation with Eternia's Design Crafted Platform

Crafted for Superior Delivery and Safety

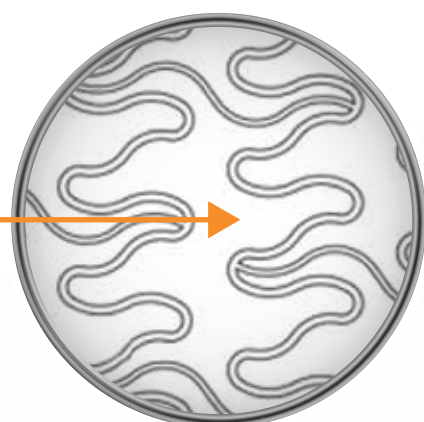

Open Cell Design

Providing enhanced side branch access, better endothelial coverage and conforms best to angulated vessels.


Proprietary WavFlo Struts

- DES with thin struts up to 60 μm
- Enhancing trackability, crossability, and radial force, while reducing the risk of restenosis and thrombosis¹



Proprietary V2V S-connectors

- Valley-to-valley long 'S' Connectors
- For higher flexibility, deliverability, Radial Strength (Force/Axial Length), and conformability¹
- Lowering stent dislodgement force to minimum 6N in 12 mm stents²

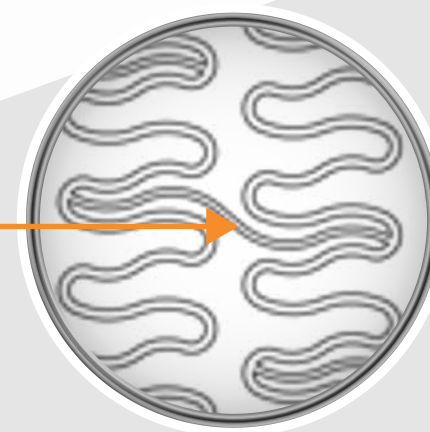

1. Kožlik M, Harpula J, Chuchra PJ, et al. Drug-Eluting Stents: Technical and Clinical Progress. Biomimetics (Basel). 2023 Feb 9;8(1):72. doi: 10.3390/biomimetics8010072.

1. Kožlik M, Harpula J, Chuchra PJ, et al. Drug-Eluting Stents: Technical and Clinical Progress. Biomimetics (Basel). 2023 Feb 9;8(1):72. doi: 10.3390/biomimetics8010072.

2. In-house data

Restore Circulation with Eternia's Precision

Specs Speak!

5-Fr Compatible Guiding Catheter

Prevents the occurrence of
bleeding, pseudoaneurysm
etc., and reduces pain at
access site

Excellent Radial Strength


Caps lesion resistance
and elastic recoil

Nearly Zero Foreshortening

Minimize operational injury
and overexpansion

Open Cell

Enhanced side
branch access

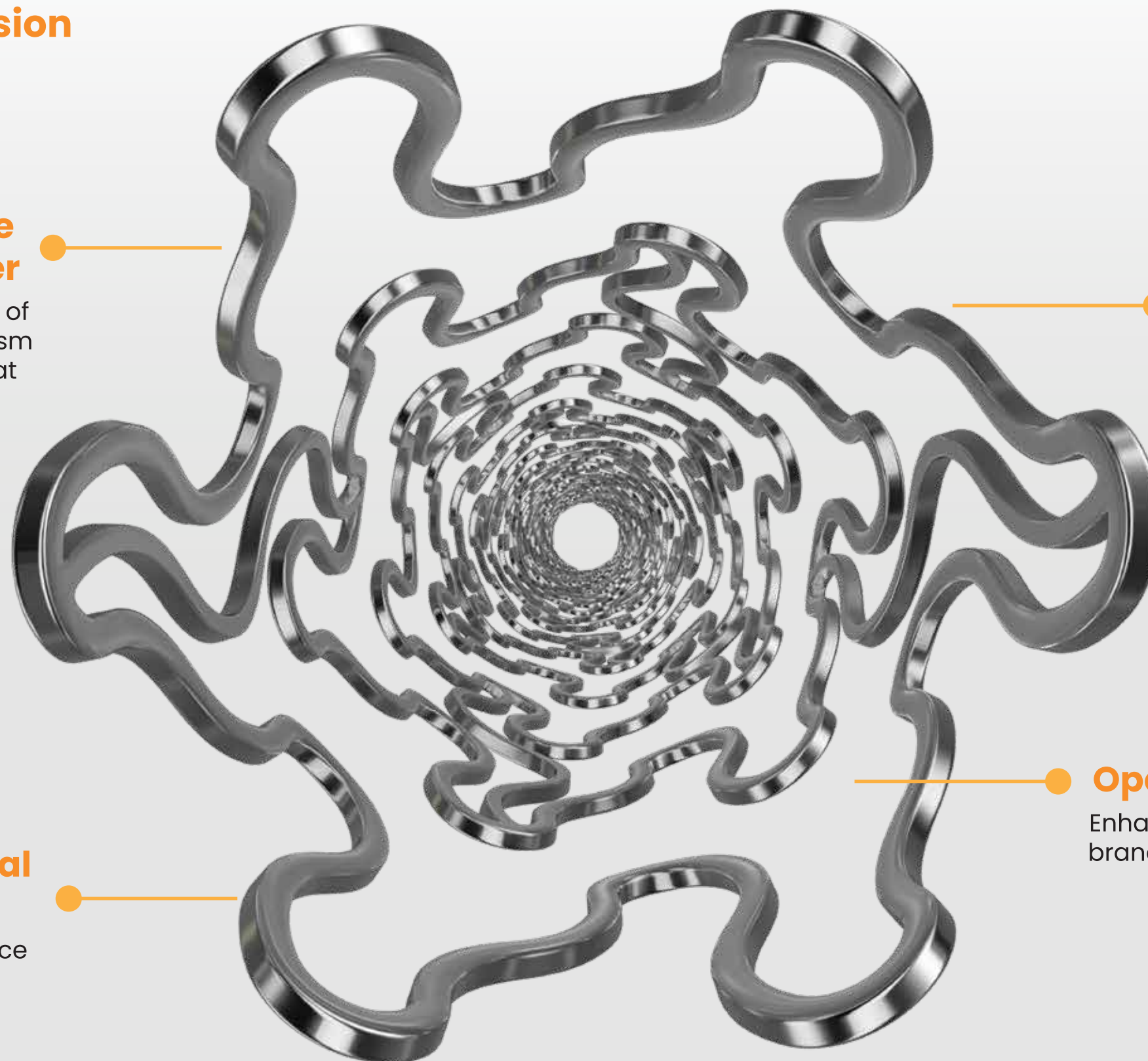

Elevate

your practice with ETERNIA, now available in a multitude of dimensions to cater to your diverse unblocking needs. Experience a new realm of possibilities for superior patient care.

ETERNIA

EVEROLIMUS ELUTING CORONARY STENT SYSTEM

Diameter (mm)	Stent Length (mm)											
	8	12	16	20	24	28	32	36	40	44	48	52
2.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.25	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.75	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
3.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
3.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
4.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
4.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓

Device Technical Specification

Design	Open Cell Design
Stent Material	L605 Cobalt Chromium
Stent Design	Open Cell Design with Valley-to-Valley S Connectors (V2V S)
Stent Strut Width	WavFlo Struts: $65 \pm 10 \mu\text{m}$ and Valley-to-Valley "S" Connectors: $55 \pm 10 \mu\text{m}$
Guide Wire Compatibility (Max)	0.014"(0.36 mm)
Polymers	Bio-degradable Polymers
Stent Strut thickness	$60 \pm 10 \mu\text{m}$
Guiding Catheter (Inner Diameter)	5 Fr Compatible
Crossing Profile	1mm ($\leq 1\text{mm}$ for 3.0 diameter)
Radial Strength (Force/Axial Length)	1.50 N/mm
Foreshortening	Nearly Zero
Flexibility	Excellent
Drug	Everolimus

Over-expansion Limit

Nominal Diameter (mm)	2.00 - 2.25 - 2.50 - 2.75 - 3.00	3.50 - 4.00 - 4.50
Number of Crowns	6	8
Post-Dilatation limit (mm)	4.00	5.50



Advant Navis Business Park Plot #7,
Office No:1012, B Wing, 10th Floor, Tower B,
Noida-Greater Noida Expy, Sector 142, Noida,
Uttar Pradesh, India 201 305.

Email: info@innvolution.com

Website: www.innvolution.com

Manufacturing License No: MFG/MD/2022/000680

Innvolution Healthcare Global Pte Ltd
Marina Bay Financial Centre Tower 1,
Level 11, 8 Marina Boulevard
Singapore 018 981.

Stent manufacturing
Plot No-105, Solitaire Industrial Park,
Village Dehmikalan, Ajmer Road,
Jaipur, Rajasthan, India 303 007.