


TECHNICAL SPECIFICATIONS

Stent length (mm)	8	12	16	20	24	28	32	36	40
Diameter (mm)	2.25	2.5	2.75	3.0	3.5	4.0			
Material	Cobalt Chromium Alloy (L-605) ASTM-F90								
Strut Thickness	73µm								
Strut Width	0.08mm (hinge) – 0.12mm (strut)								
Avg. Metal to Artery Ratio	14%								
Nominal pressure	8 ATM								
Rated burst pressure	16ATM for all sizes, except balloon diameters 4.5/5.0 (all lengths) and balloon diameters 4.0 with lengths higher than 20mm -14ATM								
Avg. foreshortening	< 10%								
Avg. recoil	< 7%								
Guiding catheter compatibility	5F for all sizes except diameters 4.5 & 5.0 - 6F								
Guiding compatibility	0.014"								

COMPLIANCE TABLE INSIGNIA

Pressure (bars)	Balloon Diameter (mm)							
	2.25	2.50	2.75	3.00	3.25	3.50	3.75	4.00
2	2.09	2.32	2.55	2.78	3.01	3.24	3.48	3.71
4	2.16	2.40	2.64	2.88	3.12	3.36	3.60	3.84
6	2.25	2.50	2.75	3.00	3.25	3.50	3.75	4.00
8	2.31	2.56	2.82	3.07	3.33	3.58	3.84	4.09
10	2.35	2.61	2.87	3.13	3.39	3.65	3.91	4.17
12	2.39	2.65	2.92	3.18	3.45	3.71	3.98	4.24
14	2.43	2.70	2.97	3.24	3.51	3.78	4.05	4.32
16	2.48	2.75	3.03	3.30	3.58	3.85	4.13	4.40
18	2.52	2.80	3.08	3.36	3.64	3.92	4.20	4.48
Nominal Pressure								
Rated Burst Pressure								

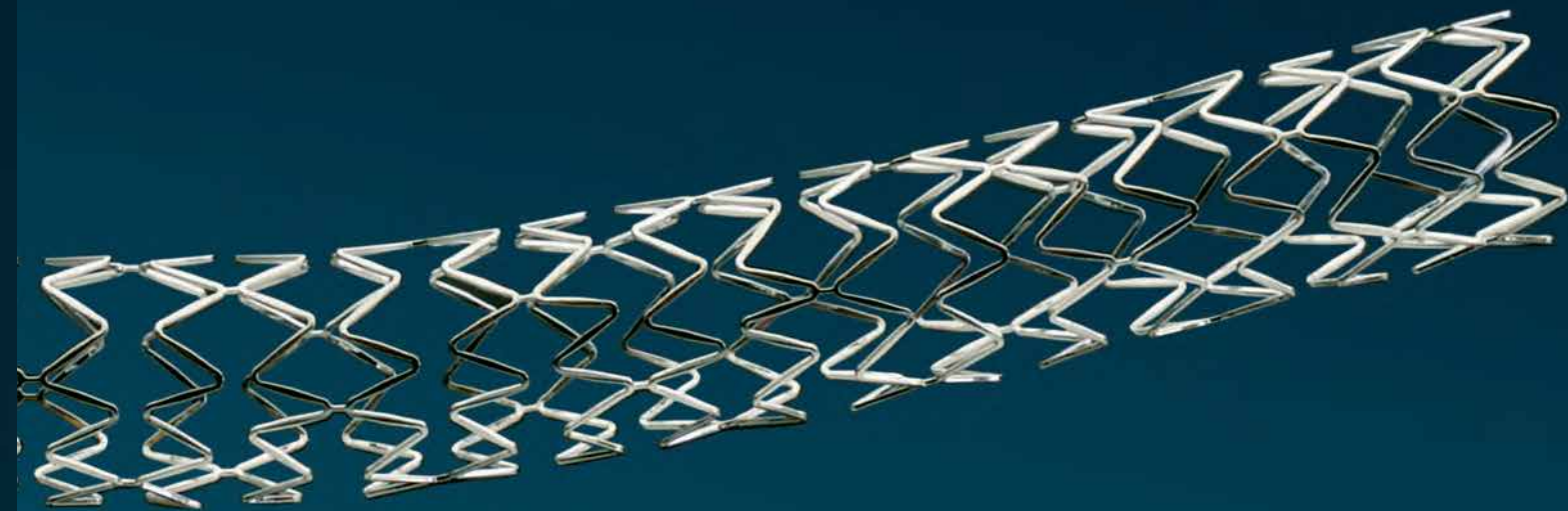

INSIGNIA - Sirolimus Drug Eluting Stent System

INSIGNIA Sirolimus Eluting Stent System elutes Sirolimus from a biodegradable polymer matrix. The drug is coated on both abluminal and luminal surfaces of the stent and is released in a controlled manner. This gives the device superior safety and efficacy.

BEST IN CLASS FEATURES

Stent & Delivery System:

- Material - L605 Cobalt Chromium
- Stent Geometry - Hybrid cell
- Low Strut Thickness- 65 Microns
- Low Balloon overhang to reduce the injury
- Low Crossing and Tip Entry Profile with Hydrophilic Coating


Drug & Polymers:

- Blend of PLA & PLGA Polymer with high Purity (>99.99%)
- Fully Biodegradable & Biocompatible Polymers
- The Lactate, produced after degradation of the polymer blend will increase the VEGF & VEGF promotes the Endothelial Cell growth
- The Loading dose of Sirolimus is 1.4 micro gram/Sq.mm surface area of the stent
- The Uniform Circumferential Coating of Drug & Polymer on the stent with a Unique Spraycoating technology maintains its integrity in in-vitro tortuous path.



Drug Release Kinetics*

Programmed Controlled Release Kinetics & Release Pattern: The initial Burst release of Sirolimus on day 1 with 20% of the loading dose inhibits the smooth-muscle cell proliferation and the rest of the Sirolimus further released from the stent inhibits both proliferation and migration of the SMC (90% of drug releases in 6 weeks).

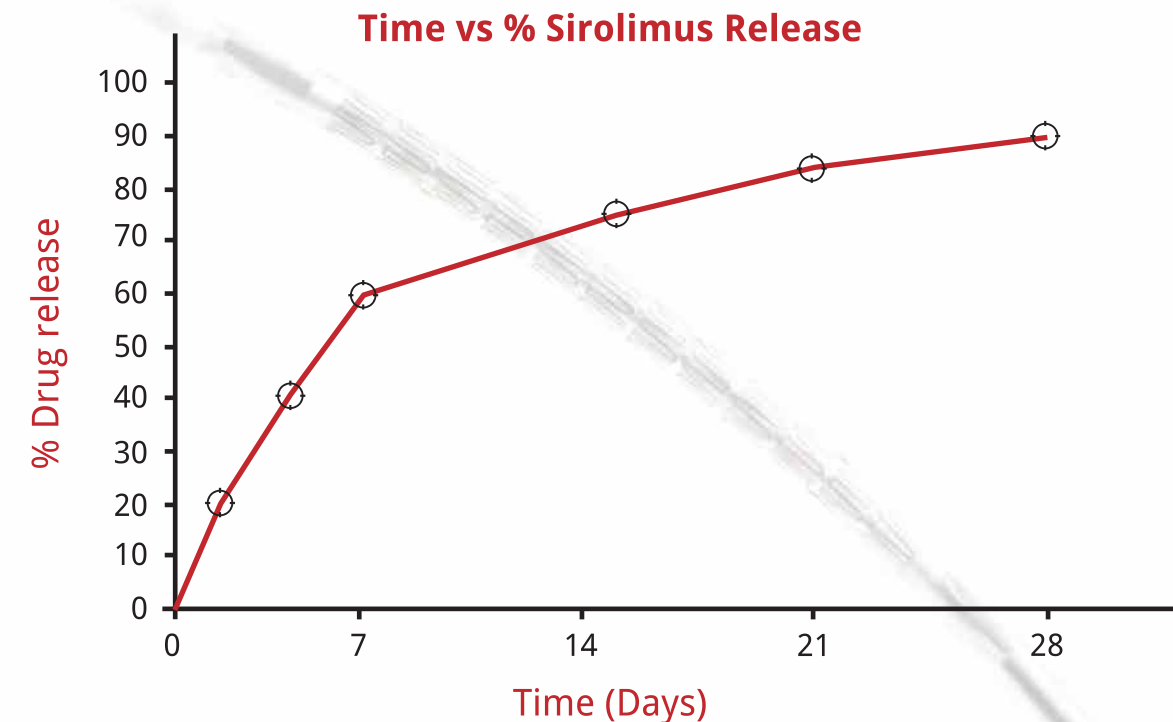

- DES to BMS: After 6 weeks there is no drug and Polymer on the surface and DES becomes BMS

Stent Apposition*

Excellent Stent Apposition as evidenced in OCT, with a score of less than 20 microns.

