

PREMIER
elite

Experience the Unmatched Flexibility of
Operating Space and Deepest Angles on a
Floor Mounted Cath Lab

CE
2265




Floor Mounted Single Plane Cath Lab

For Cardiovascular, Peripheral Vascular,
Neurovascular Diagnostic and
Interventional Procedures.






L-Arm Movement

- **$\pm 120^\circ$** Gantry movement provides unmatched space optimization enabling wide range of Cardiac, Neuro and Peripheral Vascular procedures.
- Enables all Cardiac and Neuro DSA procedures in deeper angles.
- Ease of implantation of all arrhythmia related devices.

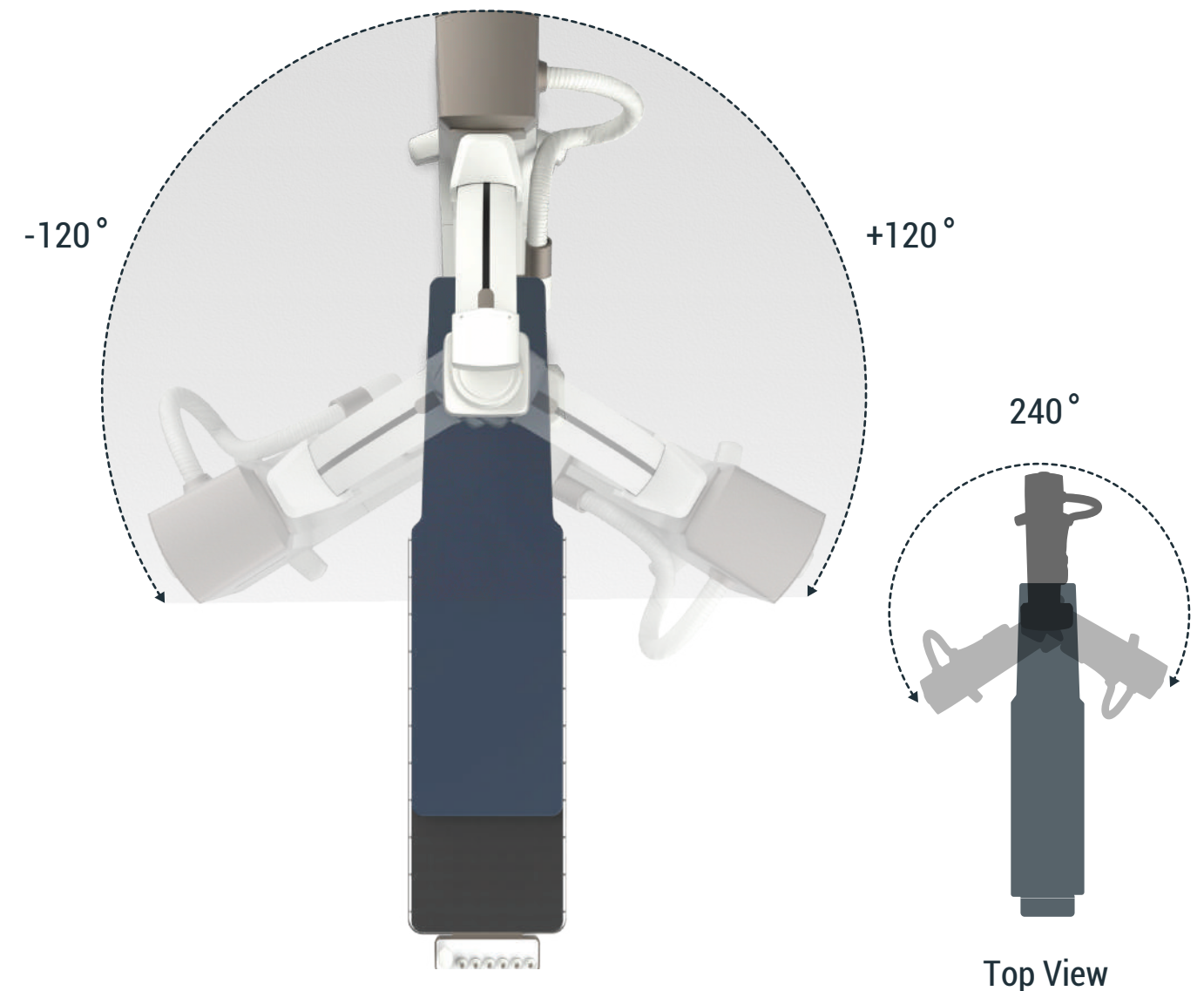

C-Arm Movement

Deeper angles ensure better access to the vasculature and enables better clinical decision making.

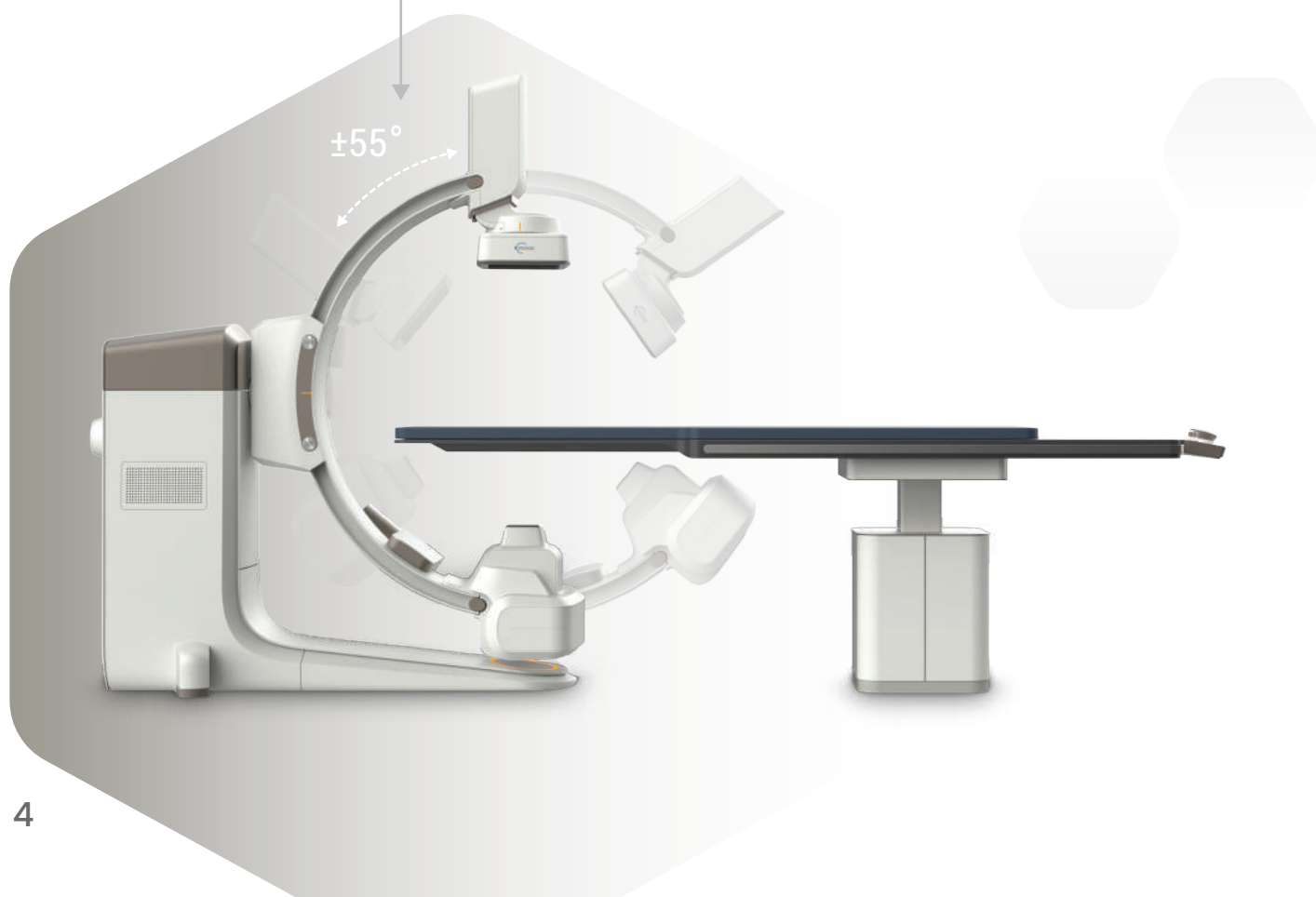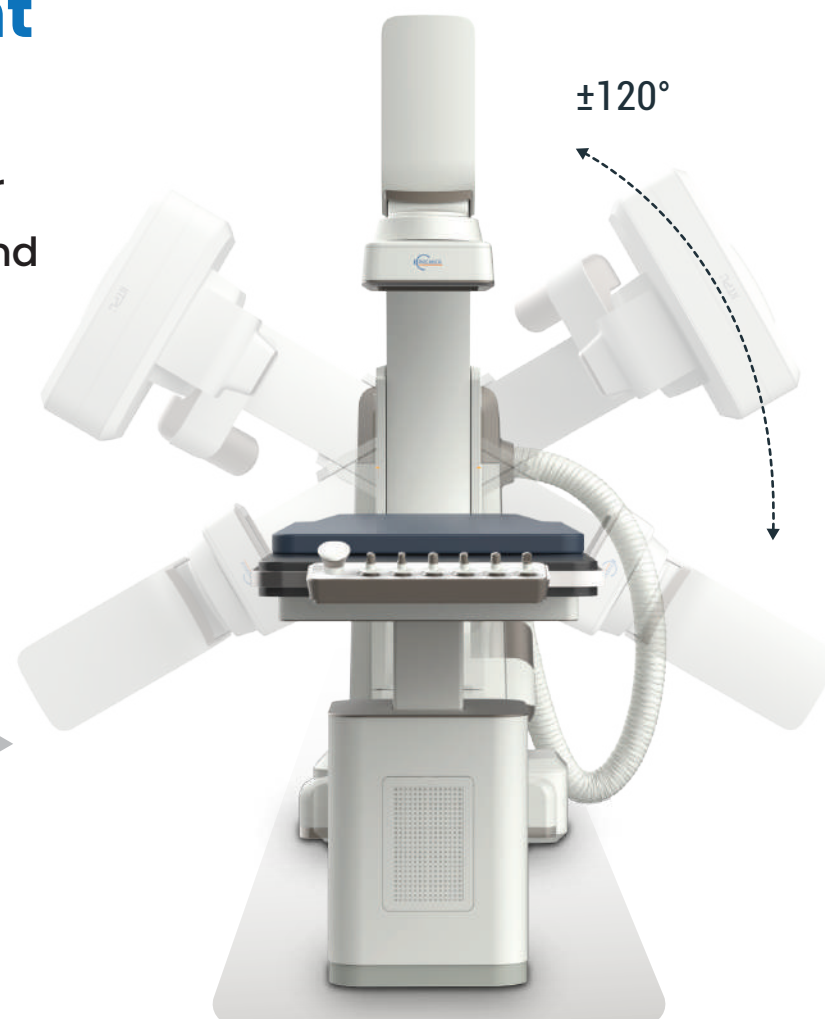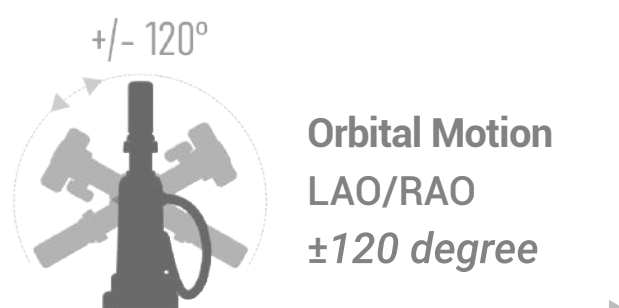

Detector Movement

- Best in class free space for the deepest angles in a floor mounted lab.
- Greater depth in detector movement enables deeper angulations thus helping better anatomical view of Heart.
- Independent Detector Movement as well as synchronized detector movement with L-Arm.

UP/Down 300 mm

Rotation ± 120 degree

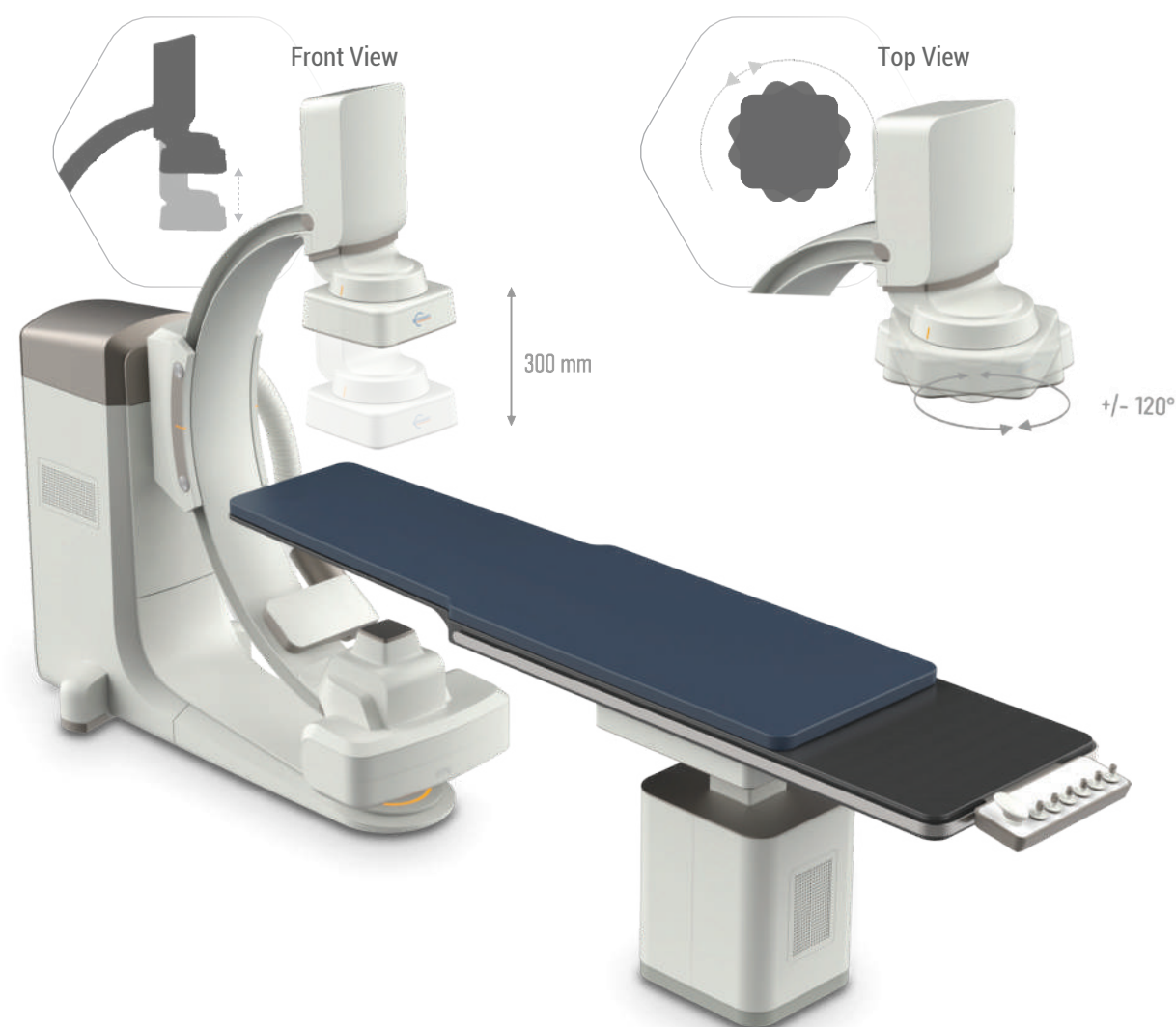


Table Swivel

Table pivot Movement of **±90 Deg**, Provides Greater flexibility to do Vascular Hybrid cases & greater access to the patient for resuscitation during emergencies.

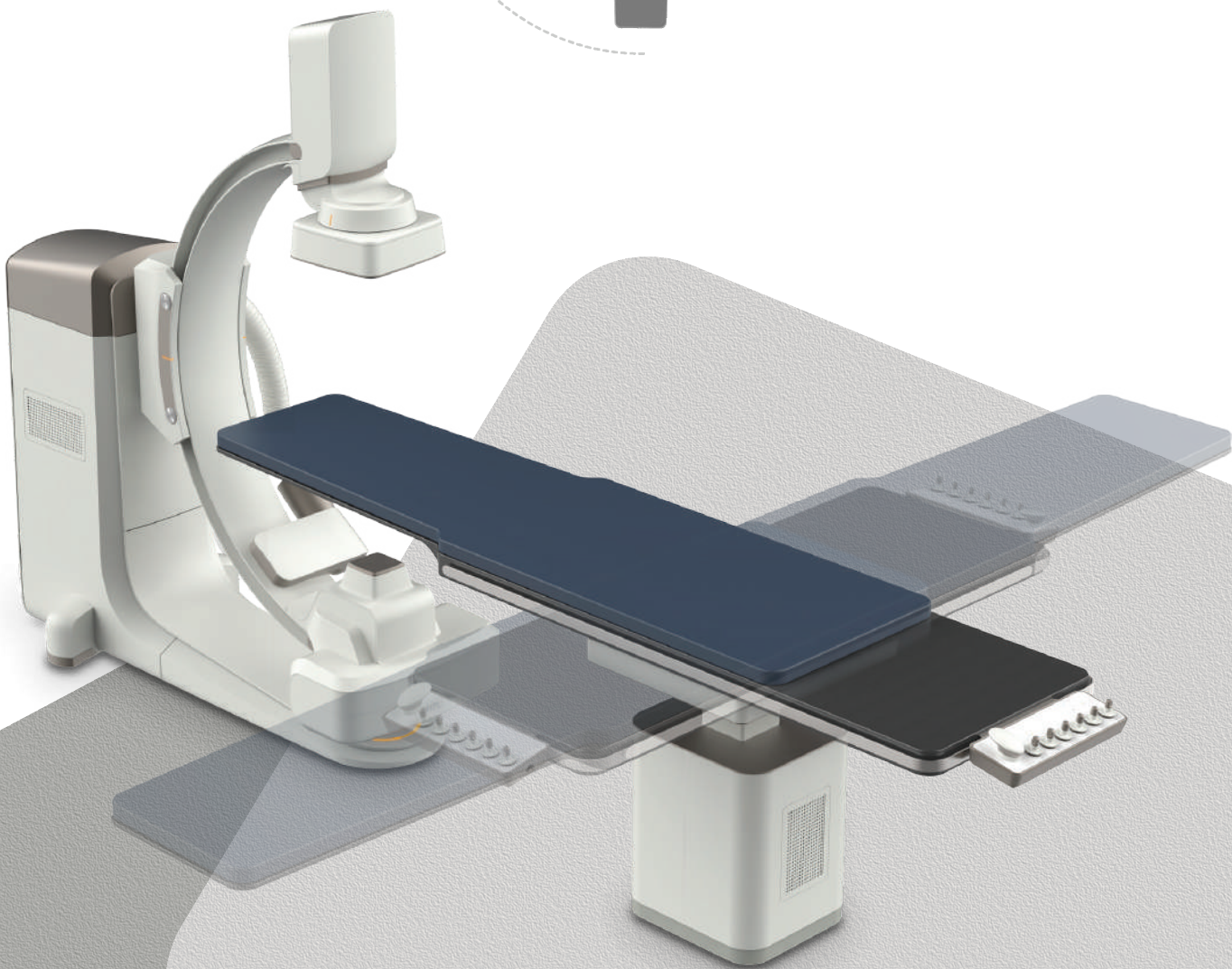
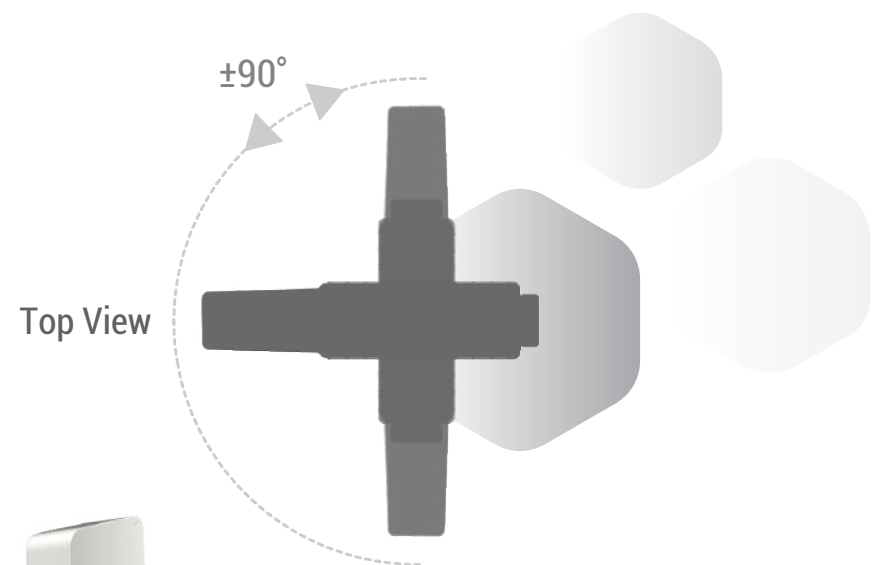


Table Sweep

Head to toe coverage to perform Peripheral Vascular, Neurovascular, Cardiac and lower extremity interventions seamlessly.


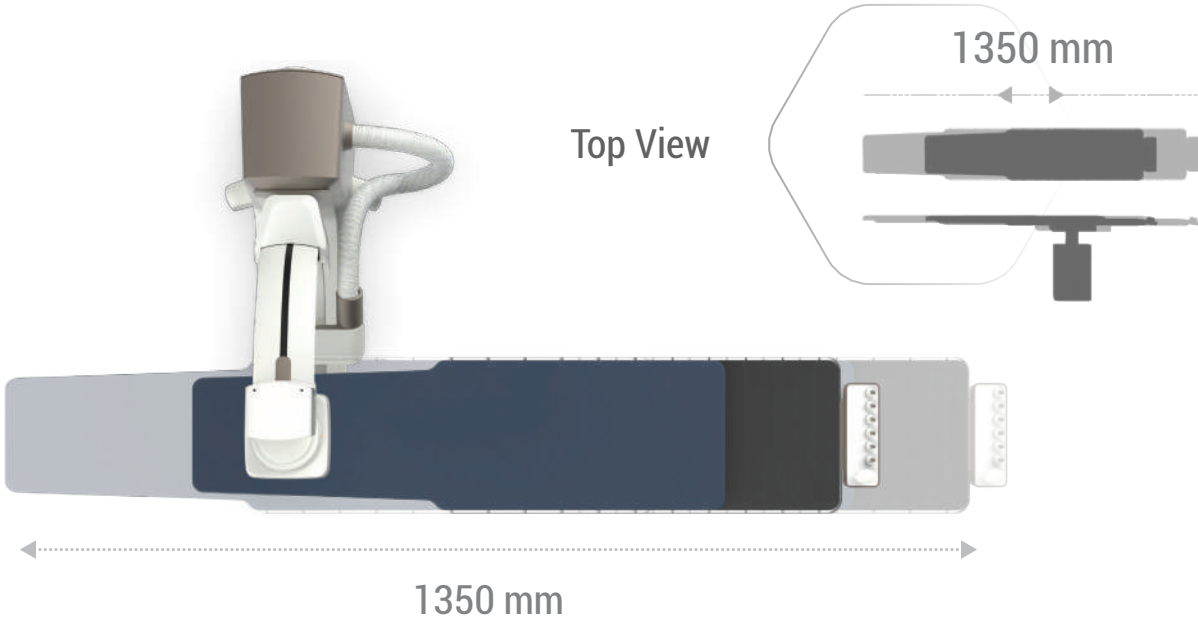



Table Height

Lower table & dismount height allows ambulatory patient to sit up without the need for a patient trolley.

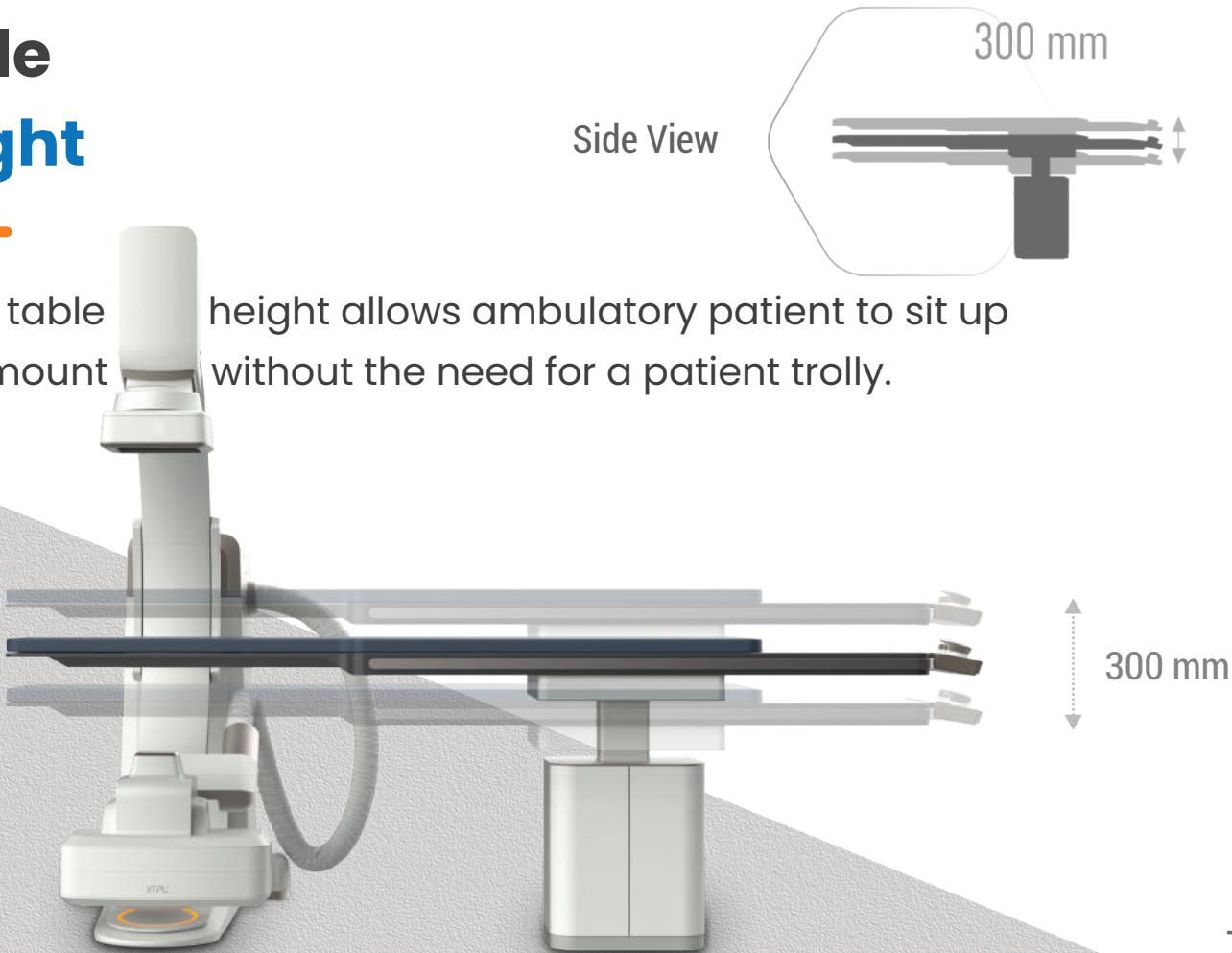

Anticollision Mechanism for Patient Safety

Non-touch based Anticollision sensors for Detector and L-Arm Side.

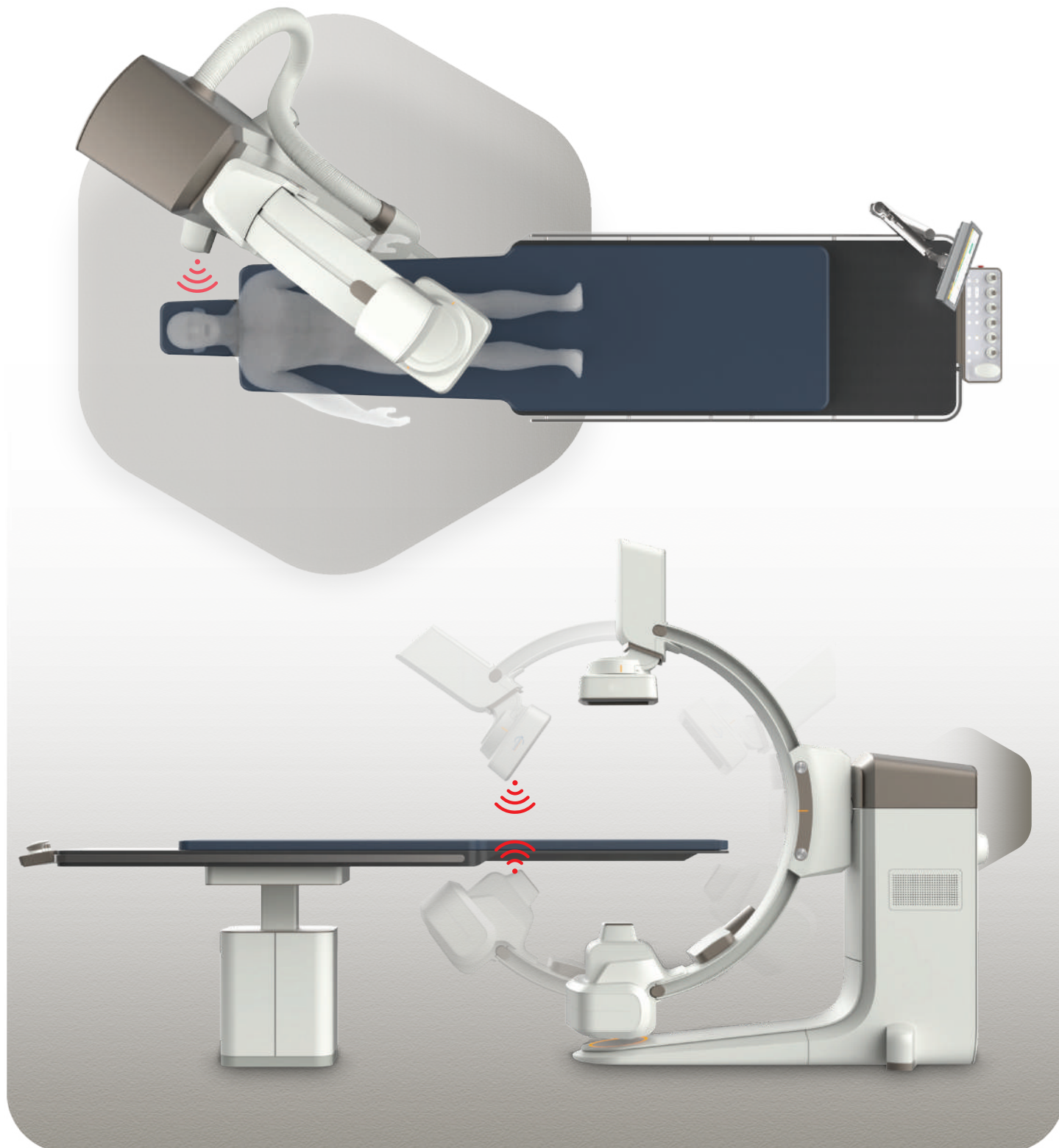

3 MHU Grid Controlled Liquid Metal Bearing Tube

Freedom to carry out lengthy and complex procedures

3 MHU High Capacity Liquid metal bearing tube supports much higher patient throughput and Continuous operation for long hours which improves tube life.



Optimized Dose Management

Grid Controlled Pulsed Fluoroscopy for sharper images and lower radiation.

Reduced Dosage

X-ray tube assembly consists of a metallic case to provide appropriate limits for X-ray.

Enhanced Image Quality


More useful radiation contributes to better image generation.

Liquid metal bearing tube is Virtually silent and delivers near real time imaging.

Smaller Focal Spots (0.4 mm & 0.8 mm) give much sharper image quality.


100 kW High Frequency X-Ray Generator

- Output power levels of **100 kW**
- **High-Level Radiography**, Continuous Fluoroscopy, Pulsed Fluoroscopy and Cineradiography.
- Better visualization of moving anatomy/OBESE Patients.
- **Higher Switching Frequency** improves image quality & reduces X-Ray dose by increasing X-Ray efficiency.
- **450 KHz** Switching Frequency enables lower high voltage ripples.

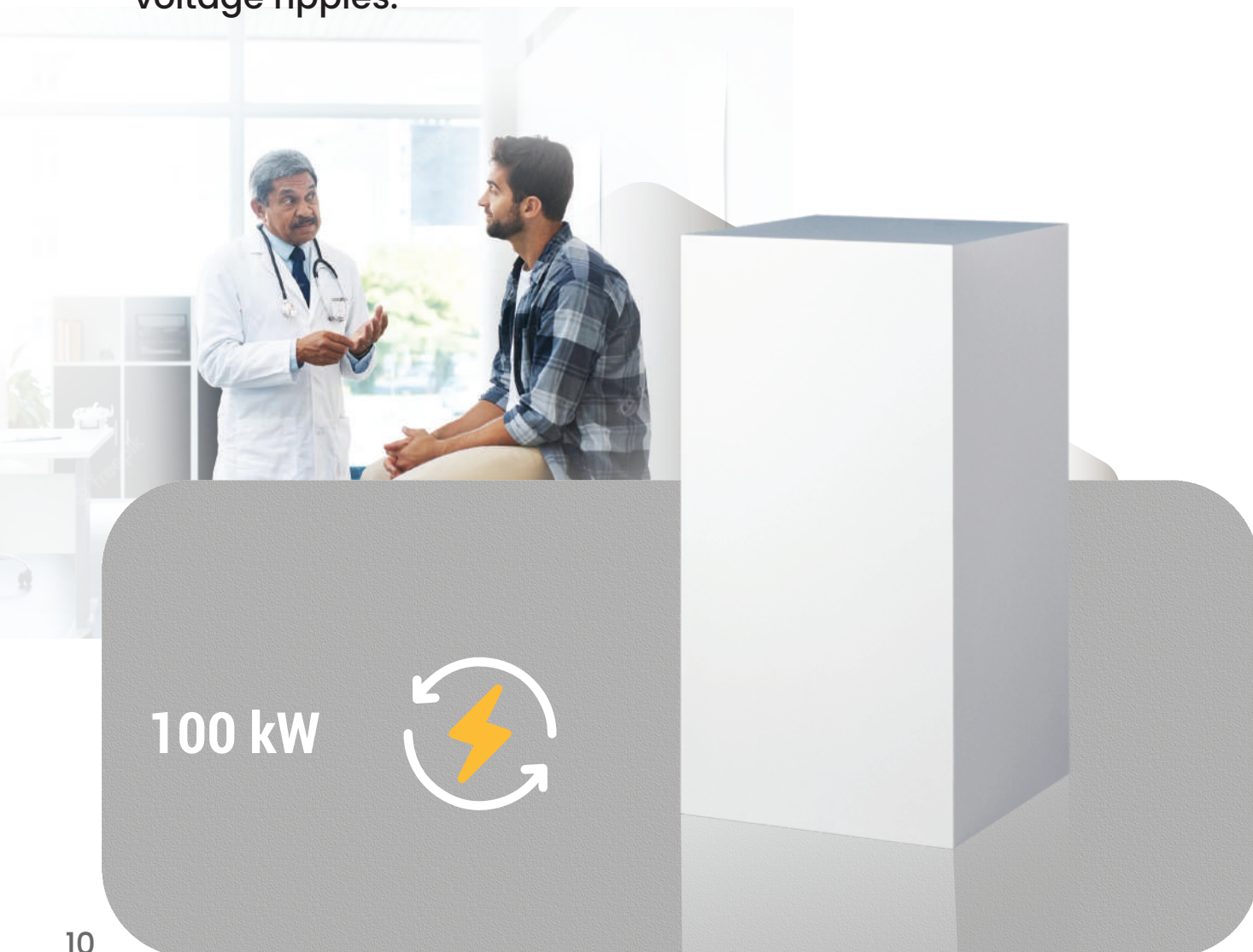

Monitor 43" inch

43" Live Monitor Display
for Sharper & Consistent
Image Quality.

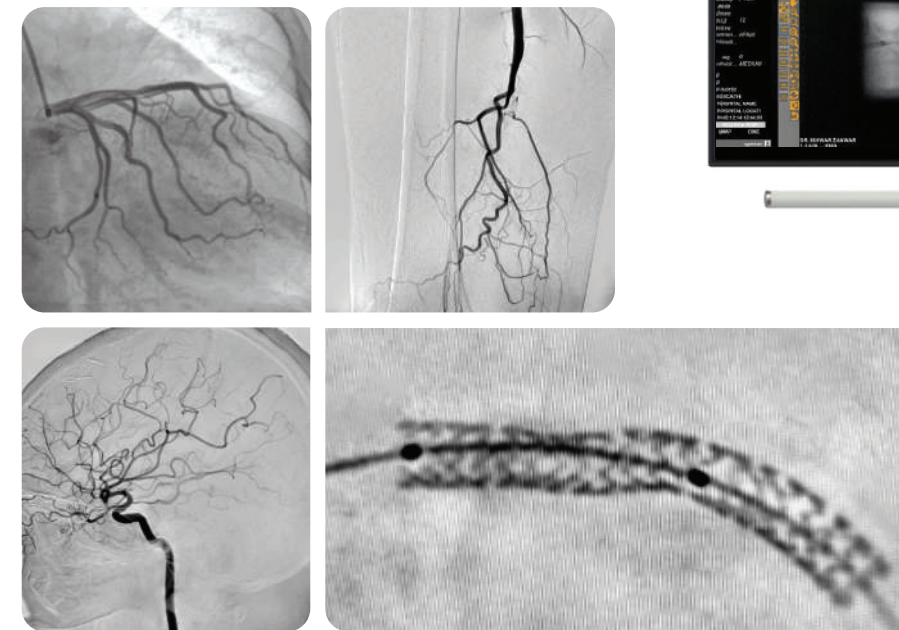

User Interface

Backlit joysticks, buttons, and
keys in the table side module for
enhanced user experience.


The table side module
includes an **ON/OFF** button
for the spot light.



Stent Clarity®


Visualise the entire stent with Stent Clarity® in high resolution.

Optional Features

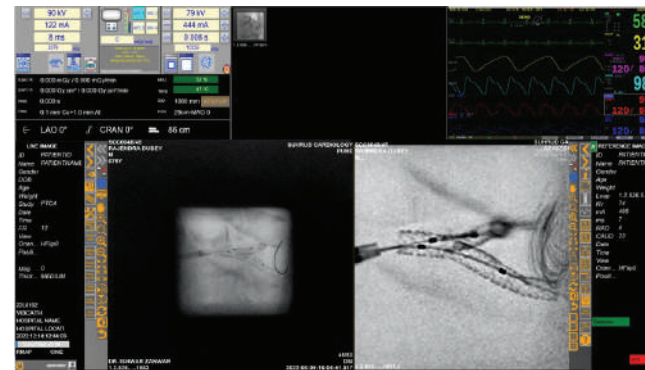

World's first AI – enabled vessel segmentation analysis, delivering results in just 2 seconds.

AI powered real-time QCA helps in automated vessel segmentation, lesion analysis, stent recommendation and automated vessel classification.

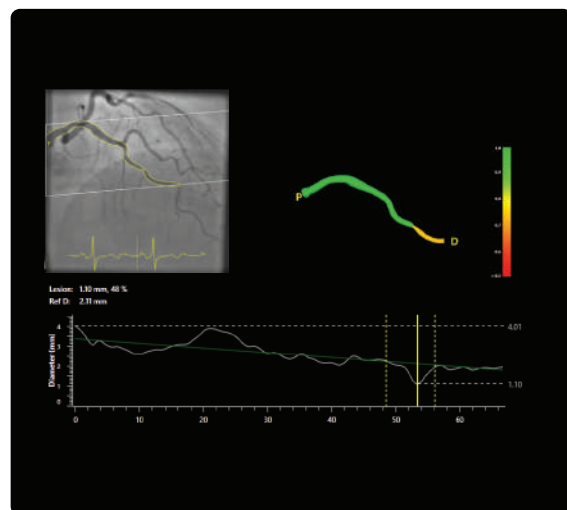

Virtual Fractional Flow Reserve

vFFR is a less invasive angio-derived FFR calculation method. Without the need of adenosine and pressure wire achieving **94%** correlation with invasive FFR measurements to enhance procedural efficiency and patient care.

Advanced Ventricular Assessment


Enhances ventricular assessment with LVA/RVA Analysis for automated workflows to precisely calculate volumes, analyze wall motion and assess ejection percentage.

Optional Features


IV-LINQ

Real-time co-registration of IVUS/OCT with angiography to obtain a detailed view of the lesion and its exact location in the coronary tree. Vendor independent Software runs with almost every X-ray, IVUS and OCT system.

IVUS Co-registration

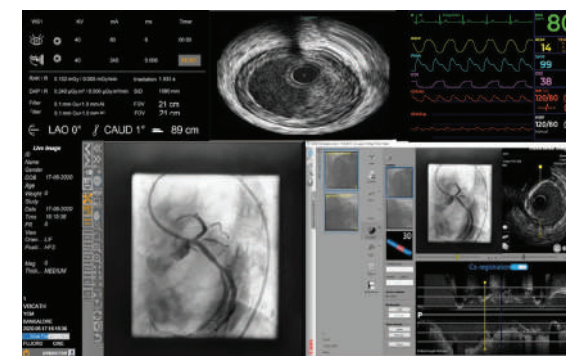

OCT Co-registration

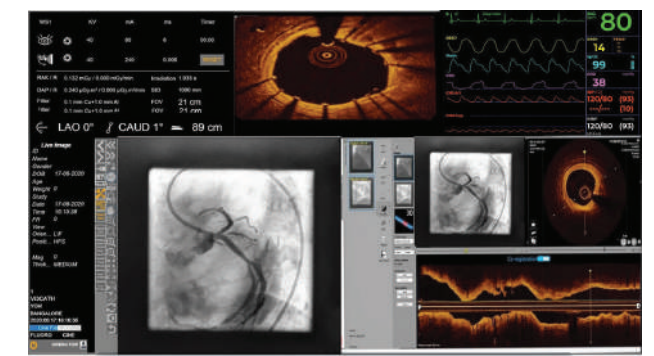

Precise Peripheral Vessel Analysis

From the abdominal aorta to carotid, renal, iliac and femoral arteries, the workflow offers unparalleled precision to improve procedural efficiency and accuracy.

Integrated Ultrasound

With combined precision and ultrasound guidance, including vascular Doppler studies, perform guided arterial punctures and transthoracic echocardiography during procedures.

Elevated TAVI Planning

(TAVI) planning before the procedure with the **3mensio Structural Heart software** that enables in elevating procedural success rates with detailed visualization and meticulous analysis.

Specification

	Gantry	
	C Arm Speed	: Up to 25 degree/Sec
	ISO Center Height	: 1110 mm
	Min. SID	: 960 mm
	Max. SID	: 1260 mm
	SID Travel	: 300 mm
	Max Free Space	: 930 mm
	L-Movement	: ±120 degree, 3 side of table
	Home Position L-arm	: Head End
	Speed	: Up to 15 degree/Sec
	2 Way Anti-collision	: Detector side and L-Arm side

	Table	
	Table top	: Carbon Fiber
	Mattress	: Memory foam
	Max. patient weight	: 200 Kg+50 kg CPR
	Min. height	: 780 mm
	Vertical stroke	: 300 mm
	Longitudinal stroke	: 1350 mm
	Table lateral stroke	: 300 mm
	Table pivot Movement	: ±90 degree

	Detector:	
	Flat Panel	: 12" diagonal
	Pixel Size	: 154 µm
	Image Matrix	: 1344x1344
	Dynamic range	: 16-bit
	DQE	: 77%

	Tube:		Generator:
	3 MHU Grid Controlled		100 kW (High frequency
	Liquid Metal Bearing tube		X-Ray generator) 40-125 kV

	Display Monitor:
	43" or 55" Live Monitor Display.
	43" 4K Industrial grade reference monitors in console room.

About us

Innvolution Healthcare is a group focused on interventional cardiology with a comprehensive portfolio of drug-eluting stents, balloon catheters, and the award-winning and fastest-growing Cath Lab Pinnacle & Premier. Innvolution Group comprises **two** companies



Innvolution Healthcare Private Limited ("IHPL"), which manufactures and distributes cardiac care consumables.



Innovation Imaging Technologies Private Limited ("IITPL"), which manufactures Cath Labs.

The company has three manufacturing facilities in India located in **Bengaluru, Vizag, and Jaipur.**

Innvolution is the fastest-growing and Asia's most awarded Cath Lab manufacturing company. The promoters of the company have a vision of creating truly differentiated solutions in Medical Imaging and offering solutions that make cardiac care more effective and accessible.



16^{*} Awards
16 National & International Awards. Asia's Most Awarded Cath Lab Company.




93%^{*}
93% customer satisfaction across 385+ installations in India.



48^{*} Hrs.
Less than 48 Hrs of downtime at any installation since inception.

* Data as of Aug 2024

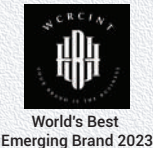



CE
2265



Address:

121F, Bommasandra Industrial Area,
Phase 1, Hulimangla Road, Hosur Main Road,
Electronic City, Bangalore 560 099, Karnataka, India

Phone:

+91 6364 882 807

Mail:

info@innvolution.com

Website:

www.innvolution.com

