



Experience Unmatched Flexibility
of Operating Space and Deepest Angles
on a Floor Mounted Cath Lab




Floor Mounted Single plane Cath Lab for Cardiovascular, Peripheral Vascular, Neuro Vascular Diagnostic and Interventional Procedures.



L-Arm Movement

± 120° Gantry movement provides unmatched Space Optimization enabling wide range of Cardiac, Neuro and Peripheral Vascular procedures.

- Enables all Cardiac and Neuro DSA procedures with deeper angles.
- Ease of implantation of all Arrhythmia related devices.

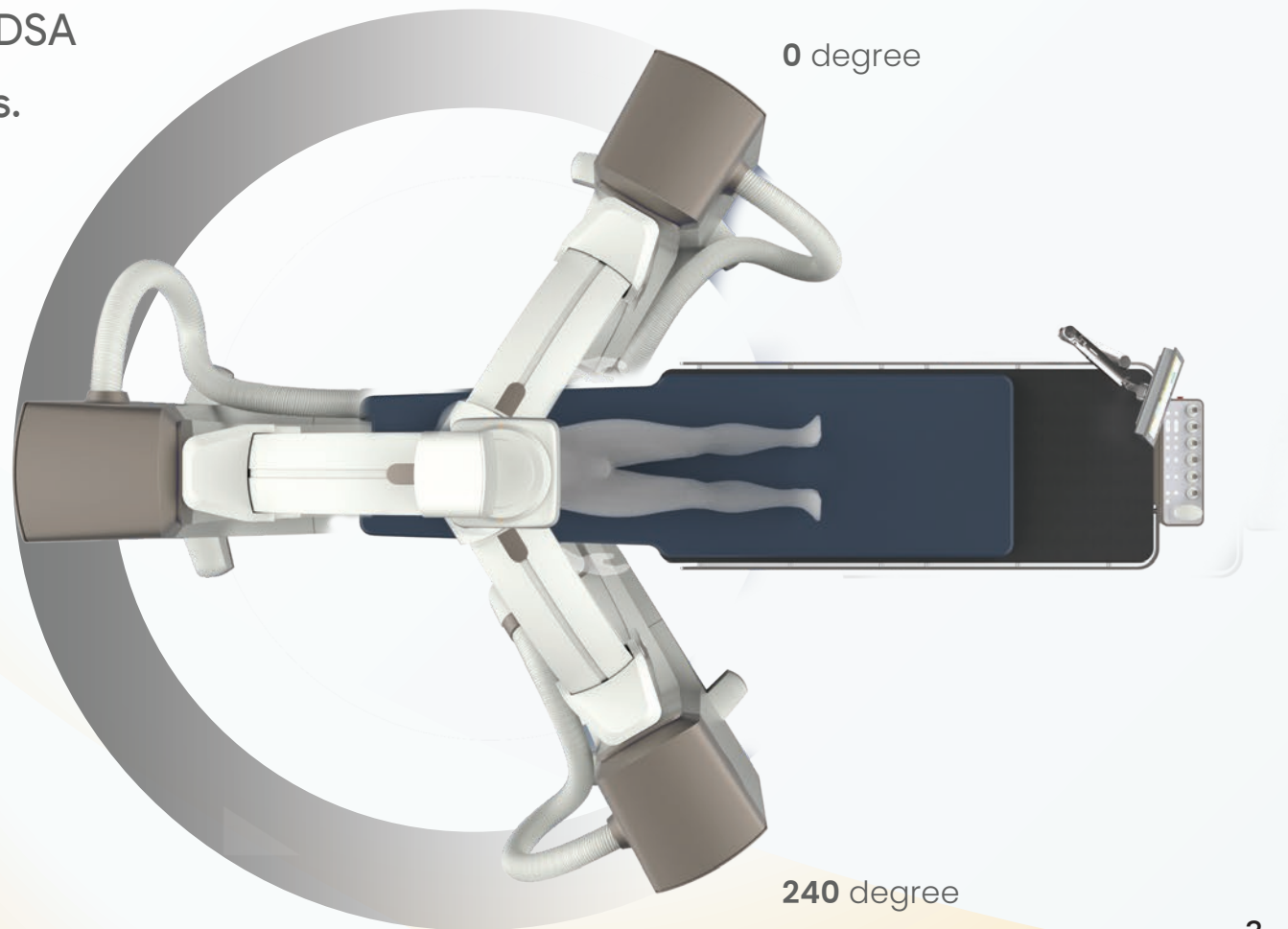

C-Arm Movement


Deeper angles ensure better access to the vasculature and enable better clinical decision making.

Angulations

CRA/CAU $\pm 55^{\circ}$



LAO/RAO $\pm 120^{\circ}$

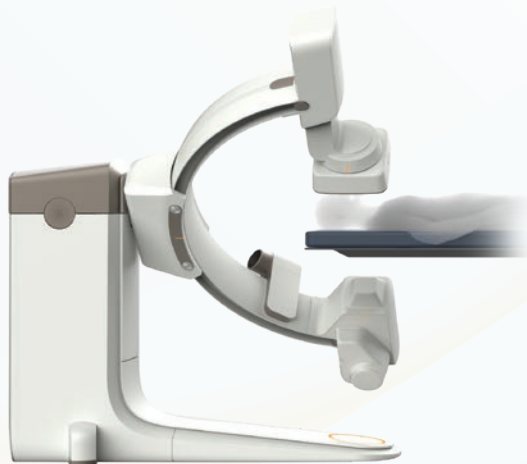

Detector Movement

- Best in class free space for the **deepest angles** in a floor mounted lab.
- Greater depth in detector movement enables deeper angulations thus helping better anatomical view of heart.
- Independent Detector Movement as well as synchronized detector movement with L-Arm.

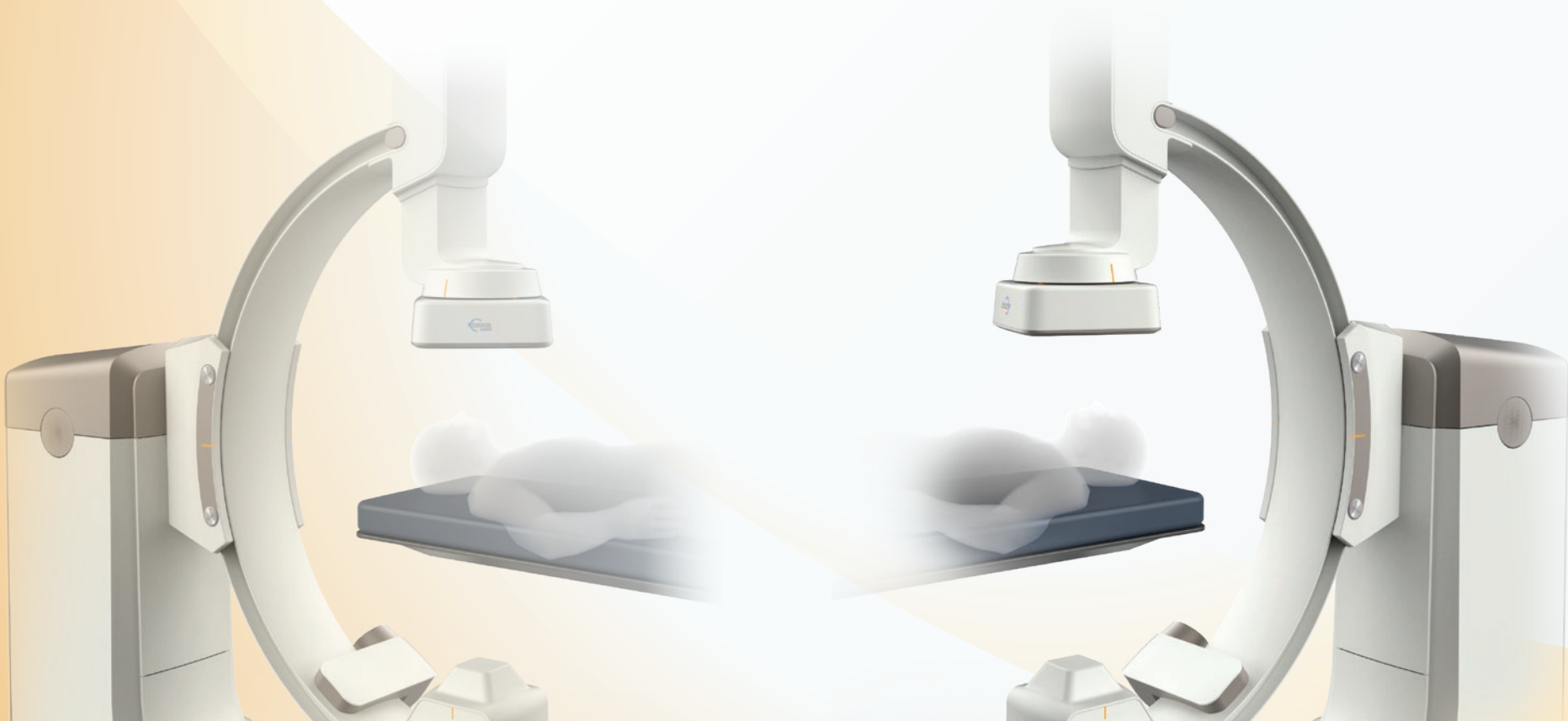


Table Swivel

Table pivot Movement of $\pm 90^\circ$,
Provides Greater flexibility to do
Vascular Hybrid cases & greater
access to the patient for resuscitation
during emergencies.


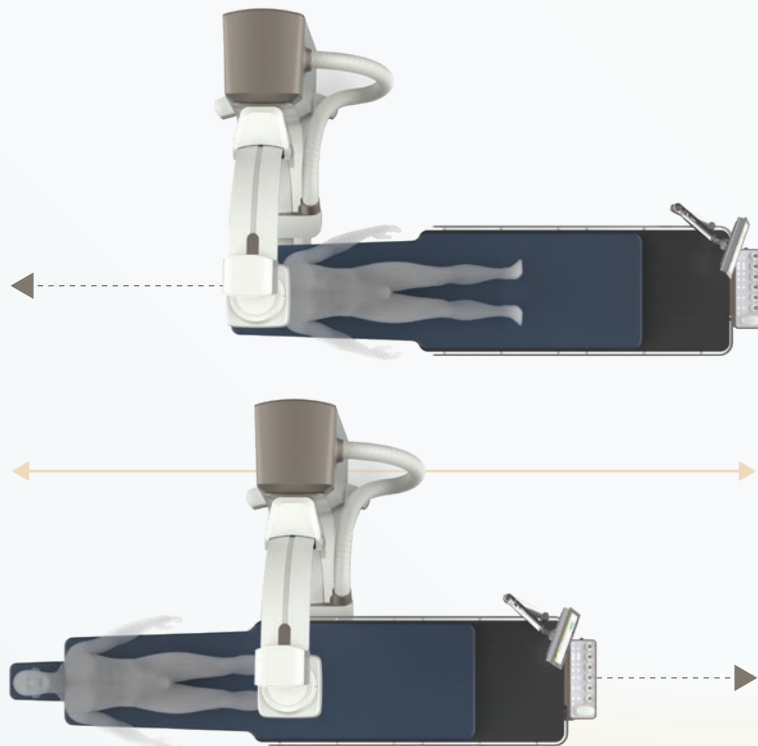


Table Sweep

Head to toe coverage to perform Peripheral
Vascular, Neurovascular, Cardiac and lower
extremity interventions seamlessly.

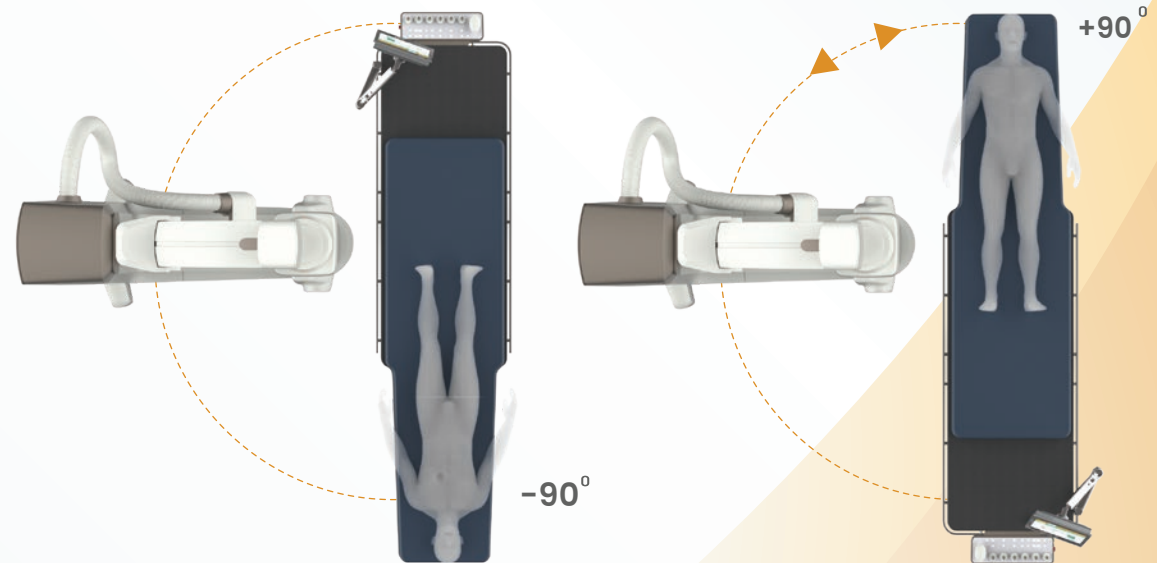

Monitor

55" Commercial Grade Monitor for Sharper and consistent Image Quality.



User Interface

Back lit Joysticks, buttons and keys in the table side module for ease of use or enhanced user experience. Spot light **ON/OFF** button provided on the table side module.




Anticollision Mechanism for Patient Safety

Non-touch based Anticollision Sensors
for Detector and Tube Side.

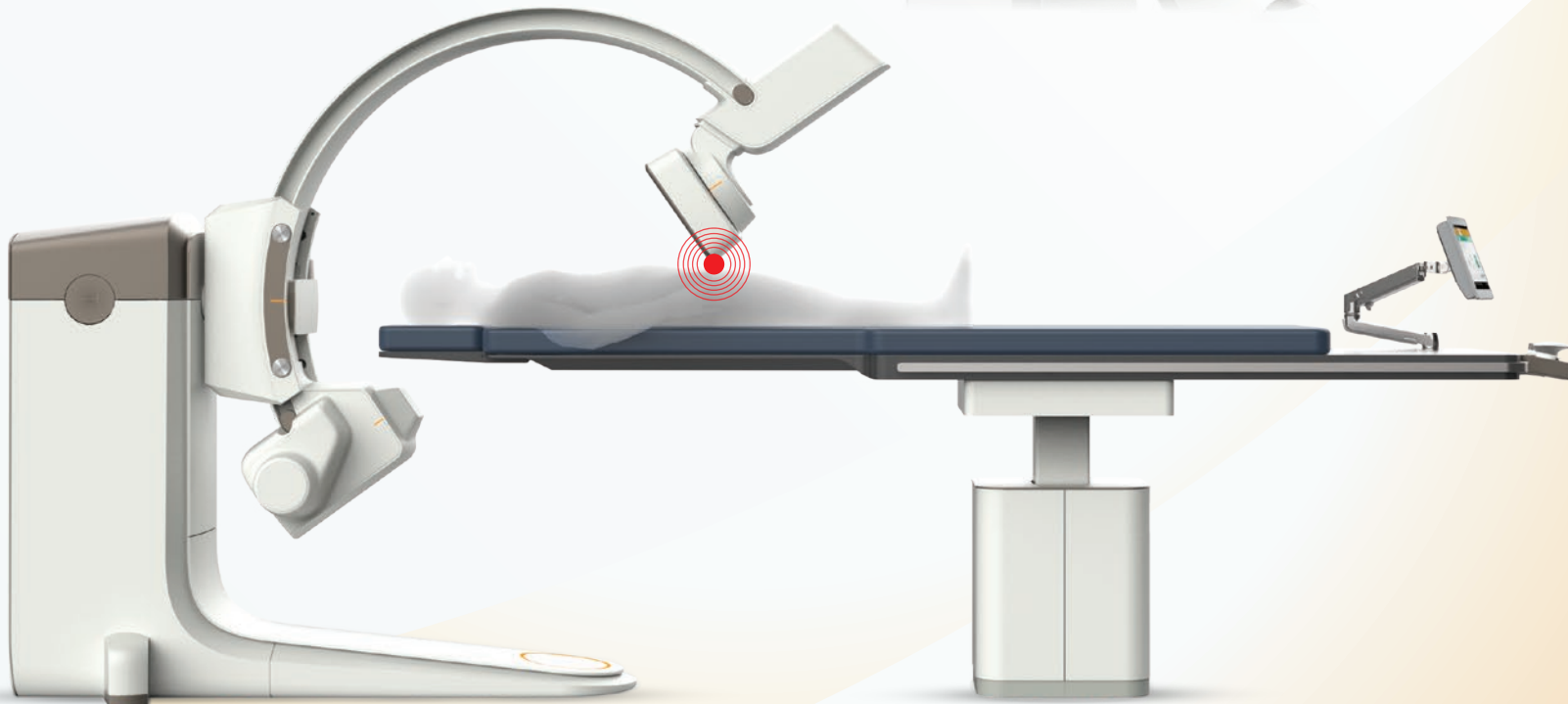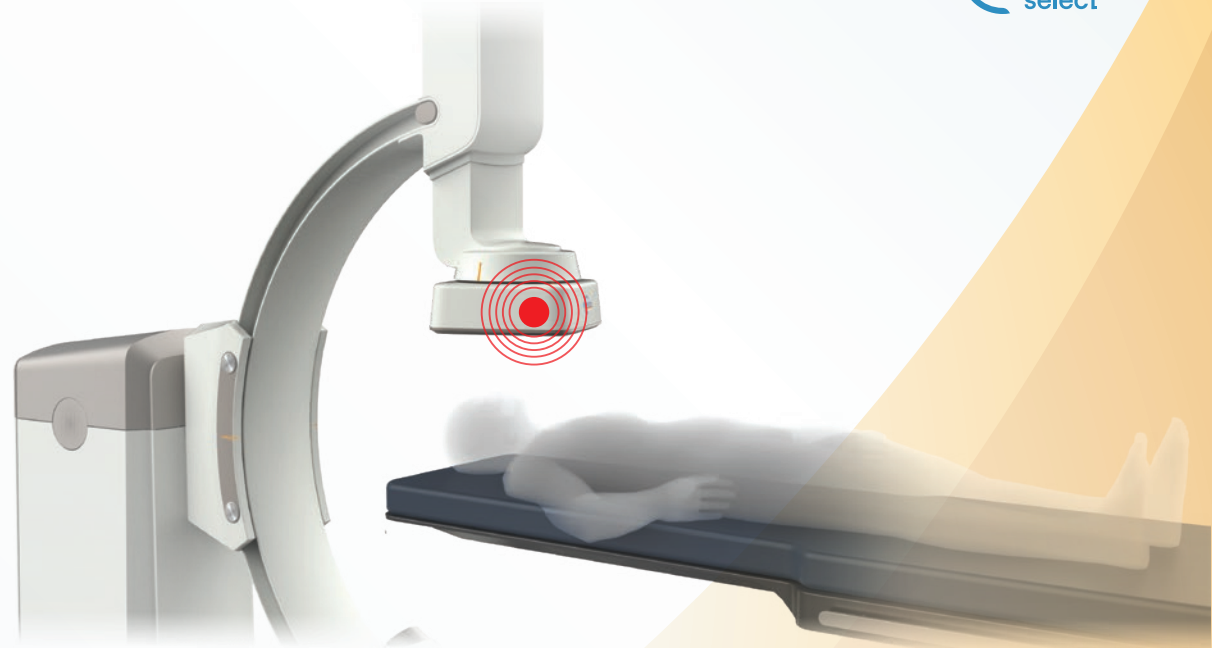

Stent Clarity®

Visualise the **entire stent** with Stent Clarity® in high resolution.

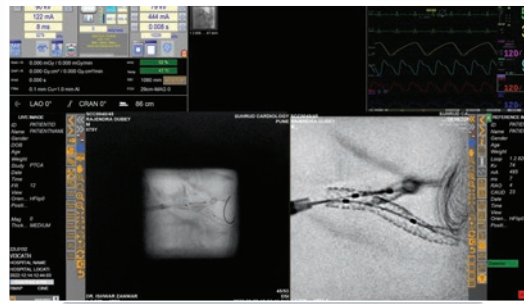


Optional Features


AINSTA QCA: World's first AI – enabled vessel segmentation analysis, delivering results in **just 2 seconds**.



AI-MAG: AI guided advanced magnification and superior resolution with **zero additional radiation exposure** and **zero additional contrast**.

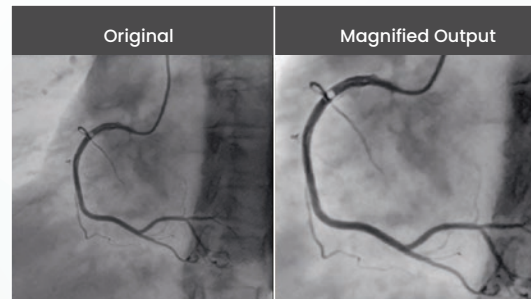

QCA: Quantitative Coronary Analysis accurately determines **vessel diameter and quantify stenosis**, which is crucial for selecting the optimal stent dimensions to suit the **original vessel diameter** and lesion length.

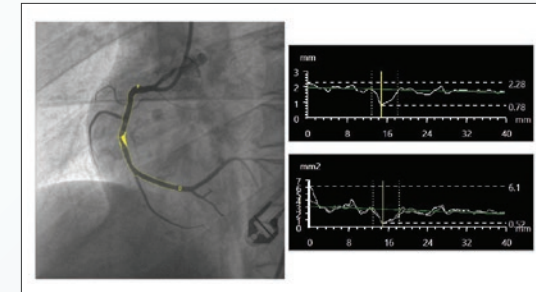

Virtual Fractional Flow Reserve

vFFR is a less invasive **angio- derived FFR calculation method**. Without the need of adenosine and pressure wire achieving **94% correlation** with invasive FFR measurements to **enhance procedural efficiency** and patient care.

IV-LINQ (Co-registration of IVUS/OCT data with angiography)

Co-register real time between angiography and IVUS/OCT with IV-LINQ to obtain **a detailed view of the lesion and its exact location in the coronary tree**.

Specification

Gantry

C-arm Speed: upto 25 degree/sec
ISO center height: 1065 mm
Min SID: 960 mm Max SID: 1260 mm
SID travel: 300 mm Max free space: 930 mm
L-arm movement: ± 120 degree
Home position of L-arm is head end
Speed: upto 15 degree/sec
2 way Anti-collision: Detector and C-arm

Tube

1.5 MHU & Oil Cooled with external Heat Exchanger
with flow and thermal safety switch.

Detector

Flat Panel 12" inch diagonal
Pixel Size 154 μ m
Image Matrix 1344x1344
16-bit Dynamic range
DQE - 77%

Table

Carbon Fiber table top
Memory foam mattress
Max. patient weight: 200 kg + 50 kg CPR
Min. height: 780 mm
Vertical stroke: 300 mm
Longitudinal stroke: 1350 mm
Table lateral stroke: 300 mm
Table pivot movement : ± 90 degree

Generator

65 kW High frequency X-Ray generator
and 40 - 125 kV

Display Monitor

55" inch Commercial Grade Display Live Monitor
One 43" inch 4K Professional grade reference monitor
and One 32" inch 4K Professional grade workstation
monitor in console room.

About us

Innvolution Healthcare is a group focused on interventional cardiology with a comprehensive portfolio of imaging and consumable products. Innvolution Group comprises two companies



Innvolution Healthcare Private Limited ("IHPL"),
which manufactures and distributes cardiac
care consumables.



Innovation Imaging Technologies Private Limited ("IITPL"),
which manufactures Cath Labs.

Innvolution's combined product portfolio includes drug-eluting stents, balloon catheters and the award-winning and fastest-growing Cath Lab Pinnacle & Premier. The company has three manufacturing facilities in India located in Bengaluru, Vizag and Jaipur.

Innvolution is the fastest-growing and Asia's most awarded Cath Lab manufacturing company. The promoters of the company have a vision of creating truly differentiated solutions in Medical Imaging and offering solutions that make cardiac care more effective and accessible.

16*

16 National and International Awards
Asia's Most Awarded Cath Lab Company

93%*

93% customer satisfaction
across **400+ installations** in India

48 Hrs*

Less than 48 Hrs of downtime
at any installation since inception

* Data as of Oct 2024

