


Flex your Way into Complex
Lesions with

SIROFLEX

SIROLIMUS ELUTING CORONARY STENT SYSTEM

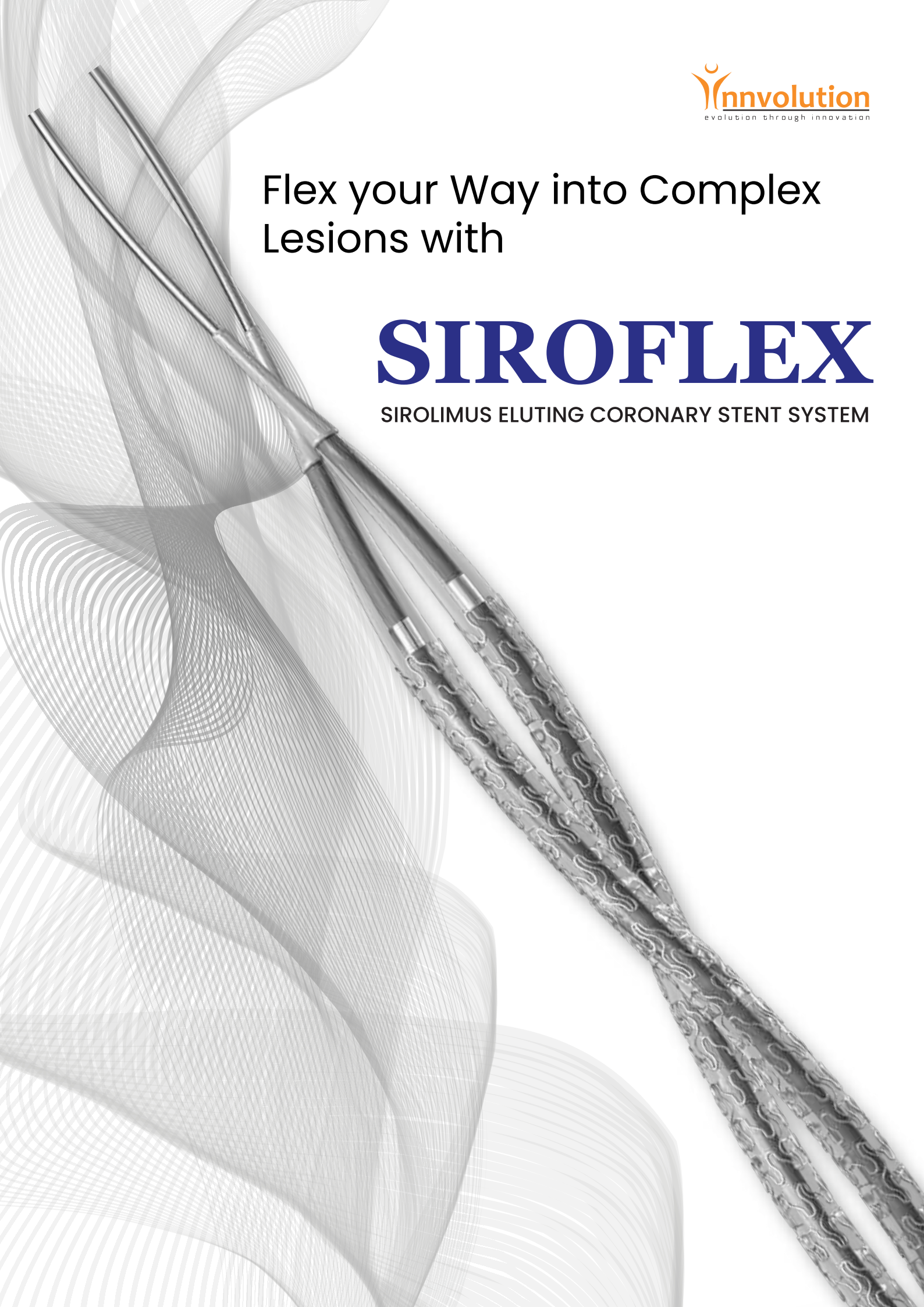





An Innovative Sirolimus-Eluting Coronary Stent System Extreme Flexibility for

Proprietary

- Valley-to-valley long 'S'Connectors
- For higher Flexibility, Deliverability, Radial Strength (Force/Axial Length) and Conformability
- Lowering stent dislodgement force to minimum 6N in 12 mm stents

Proprietary WavFlo Struts


- DES with thin struts up to 60 μm
- Enhancing trackability, crossability, and radial force, while reducing the risk of restenosis and thrombosis

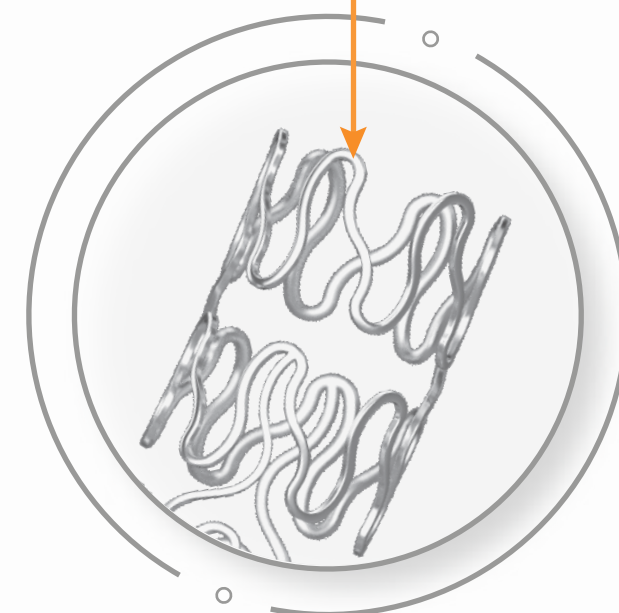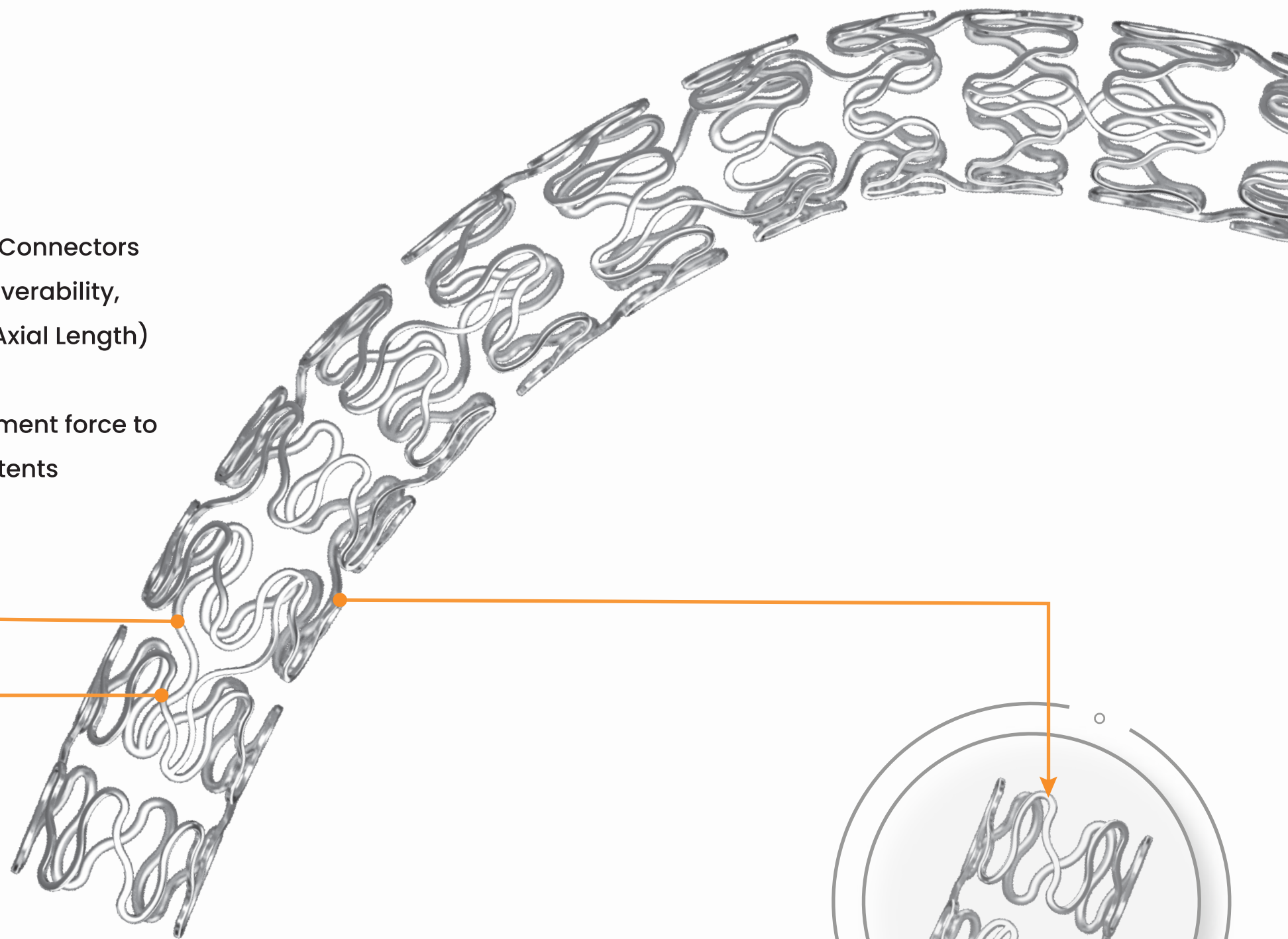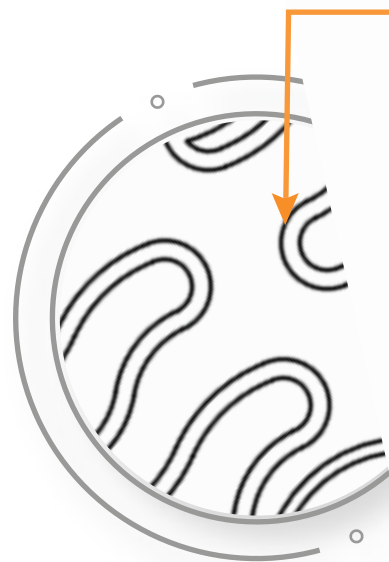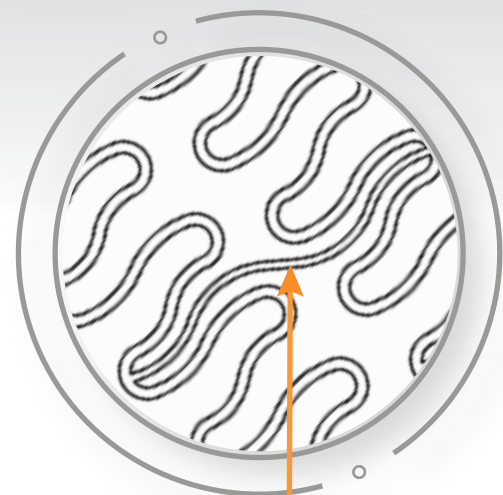

An Innovative Sirolimus Stent with Extreme Flexibility for Complex Lesions


Spacious open cell design

- Helps navigate tortuous and serpentine vessels
- Providing high flexibility
- Providing enhanced side branch access

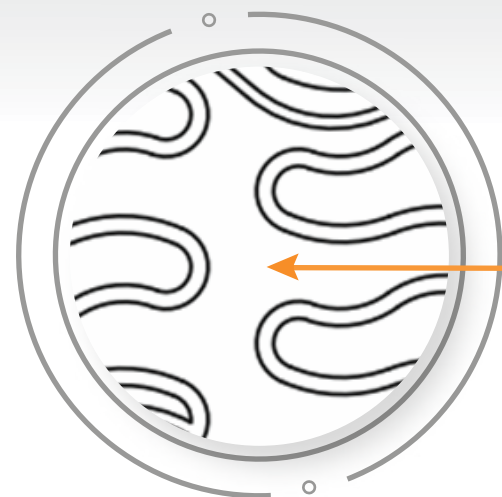



Enhanced deliverability with 60 µm ultrathin struts

- Preventing stent thrombosis & restenosis
- Reducing deep wall trauma

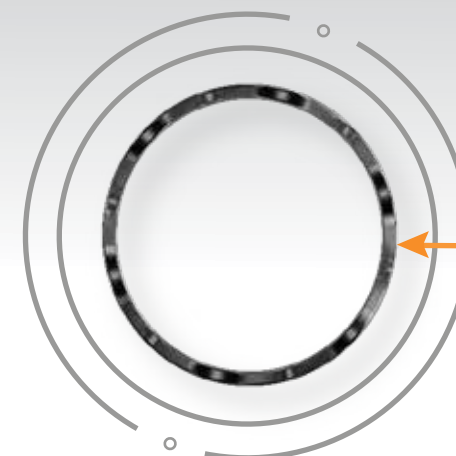


Perfectly engineered smooth edges

- Preventing artery injury or dissection

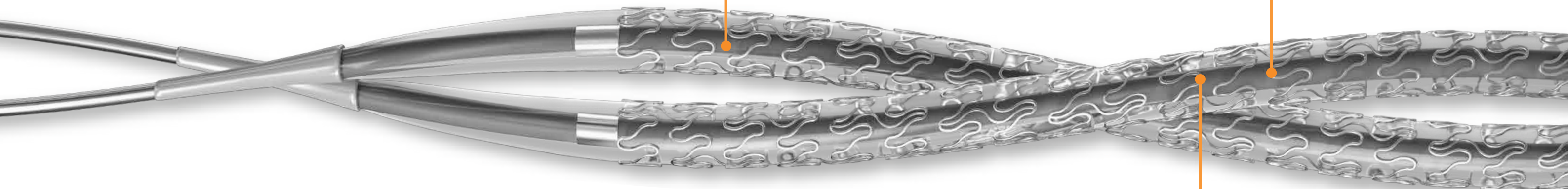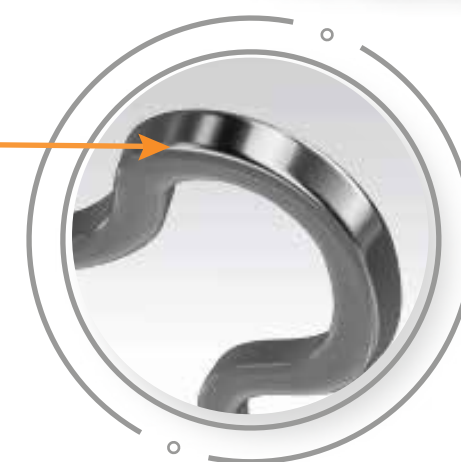


With Technology Translating to Effective Clinical Outcomes

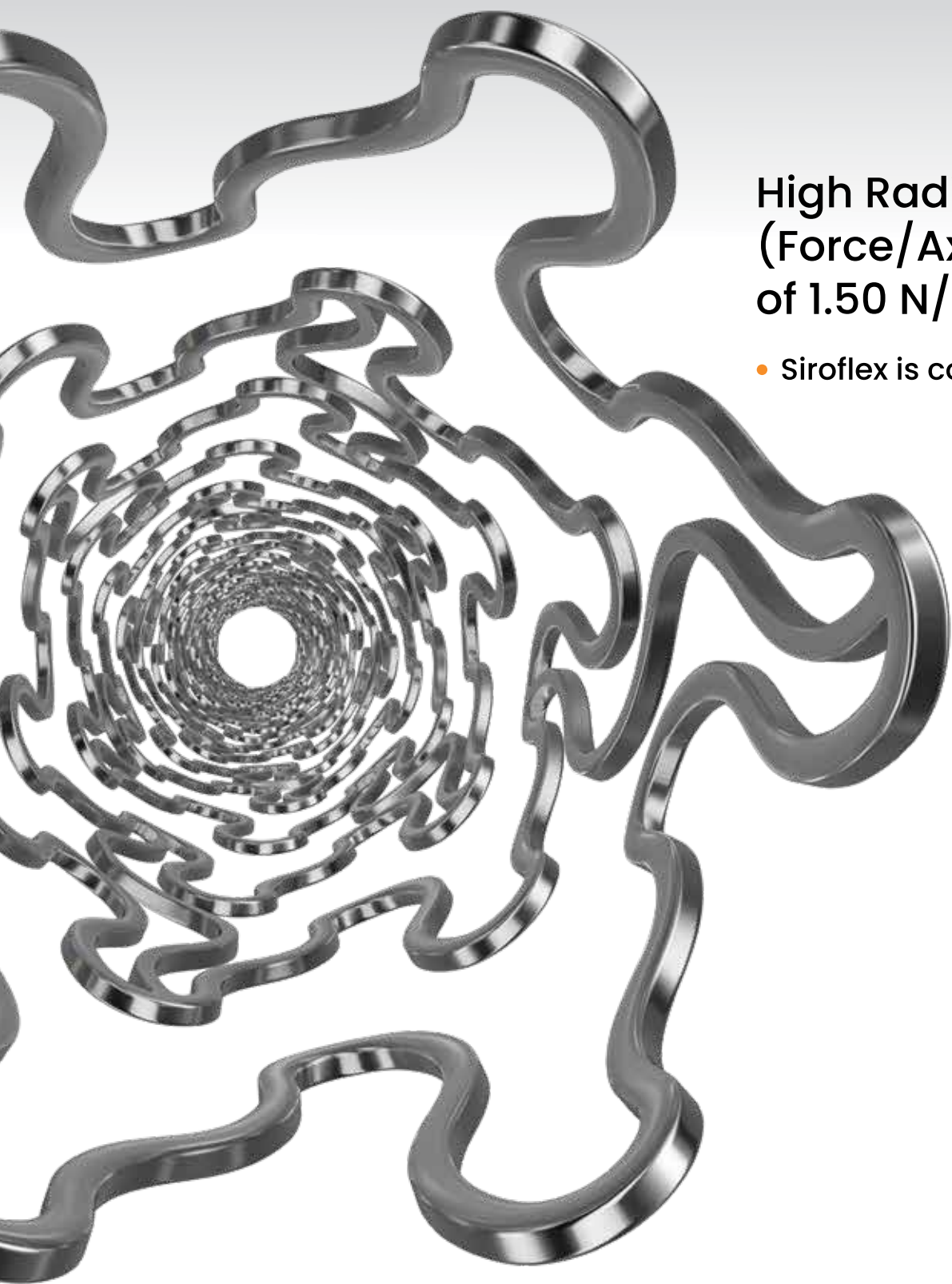

**High Radial Strength
(Force/Axial Length)
of 1.50 N/mm**

- Siroflex is collapse resistant

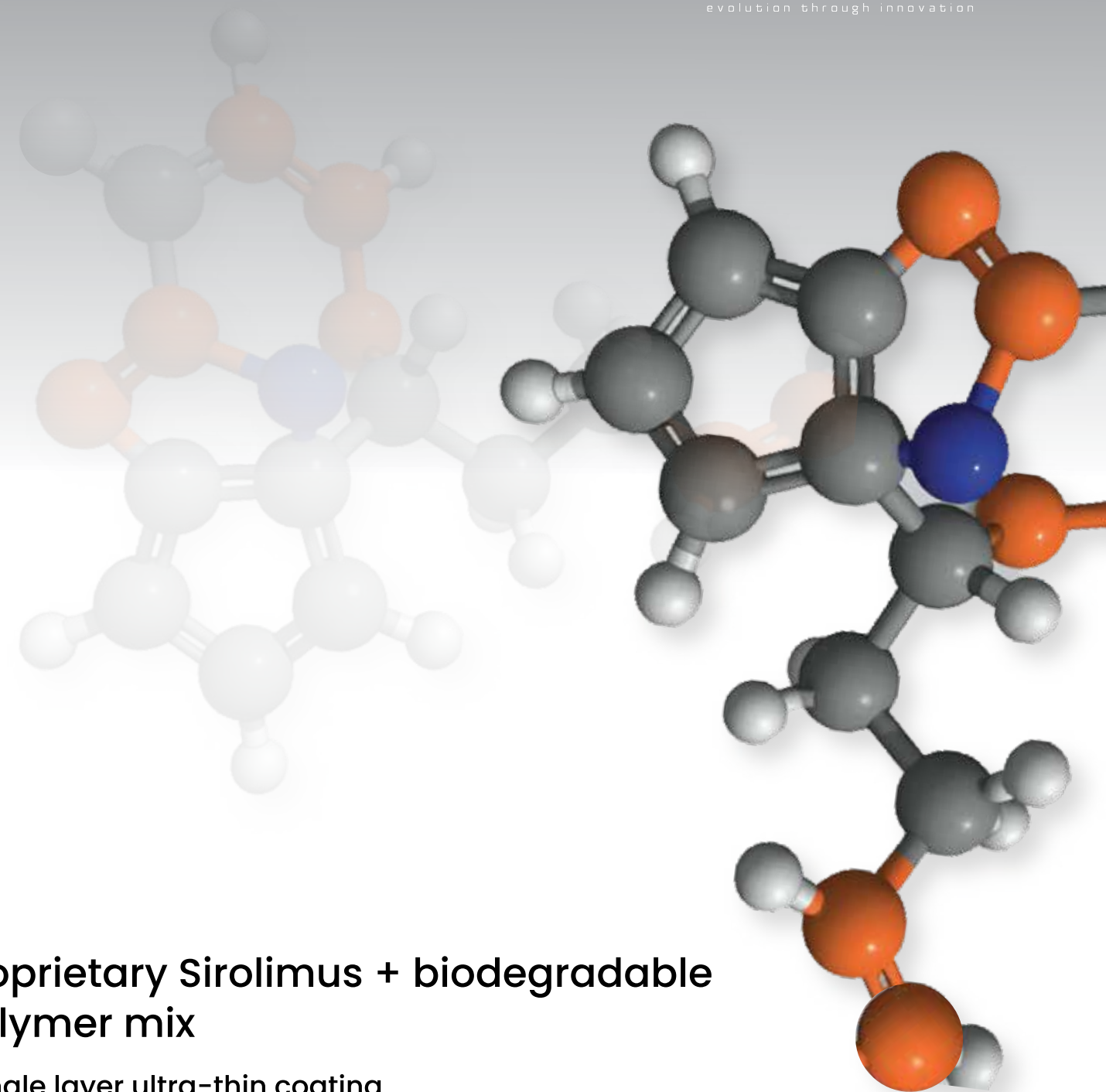

Proprietary Sirolimus + biodegradable polymer mix

- Single layer ultra-thin coating
- Decreases tissue-polymer contact
- Reduces late-stent thrombosis

SIROFLEX

SIROLIMUS ELUTING CORONARY STENT SYSTEM

Navigate through complex lesions effortlessly with Siroflex. Offering a range of dimensions to accommodate your diverse patient requirements.

Diameter (mm)	Stent Length (mm)											
	8	12	16	20	24	28	32	36	40	44	48	52
2.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.25	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
2.75	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
3.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
3.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
4.00	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
4.50	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓

Device Technical Specification

Design	Open Cell Design
Stent Material	L605 Cobalt Chromium
Stent Design	Open Cell Design with Valley-to-Valley S Connectors (V2V S)
Stent Strut Width	WavFlo Struts: $65 \pm 10 \mu\text{m}$ and Valley-to-Valley "S" Connectors: $55 \pm 10 \mu\text{m}$
Guide Wire Compatibility (Max)	0.014"(0.36 mm)
Polymers	Bio-degradable Polymers
Stent Strut thickness	$60 \pm 10 \mu\text{m}$
Guiding Catheter (Inner Diameter)	5 Fr Compatible
Crossing Profile	1mm ($\leq 1\text{mm}$ for 3.0 diameter)
Radial Strength (Force/Axial Length)	1.50 N/mm
Foreshortening	Nearly Zero
Flexibility	Excellent
Drug	Sirolimus

Over-expansion Limit

Nominal Diameter (mm)	2.00 - 2.25 - 2.50 - 2.75 - 3.00	3.50 - 4.00 - 4.50
Number of Crowns	6	8
Post-Dilatation limit (mm)	4.00	5.50

SiroFlex: SFX-BRC-VI-02012024



Advant Navis Business Park Plot #7,
Office No:1012, B Wing, 10th Floor, Tower B,
Noida-Greater Noida Expy, Sector 142, Noida,
Uttar Pradesh, India 201 305.

Email: info@innvolution.com

Website: www.innnvolution.com

Manufacturing License No: MFG/MD/2022/000680

Innnvolution Healthcare Global Pte Ltd
Marina Bay Financial Centre Tower 1,
Level 11, 8 Marina Boulevard
Singapore 018 981.

Stent manufacturing
Plot No-105, Solitaire Industrial Park,
Village Dehmikalan, Ajmer Road,
Jaipur, Rajasthan, India 303 007.